


**Terfenadine****Column:** N/A**Dimensions:** N/A**Order No:** N/A**Elution Type:** N/A**Eluent A:** Hexane/Dichloroethane/Ethanol-TFA (58:35:7) w/ ethanol-TFA (20:1)**Gradient**      **Step No.**      **Time (min)**      **Pct A****Profile:**              **1**              0              100**Flow Rate:** 1 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 260 nm (ambient)**ANALYTES:****1** Terfenadine

Retention Time: 11.21 min

**2** Terfenadine

Retention Time: 11.58 min

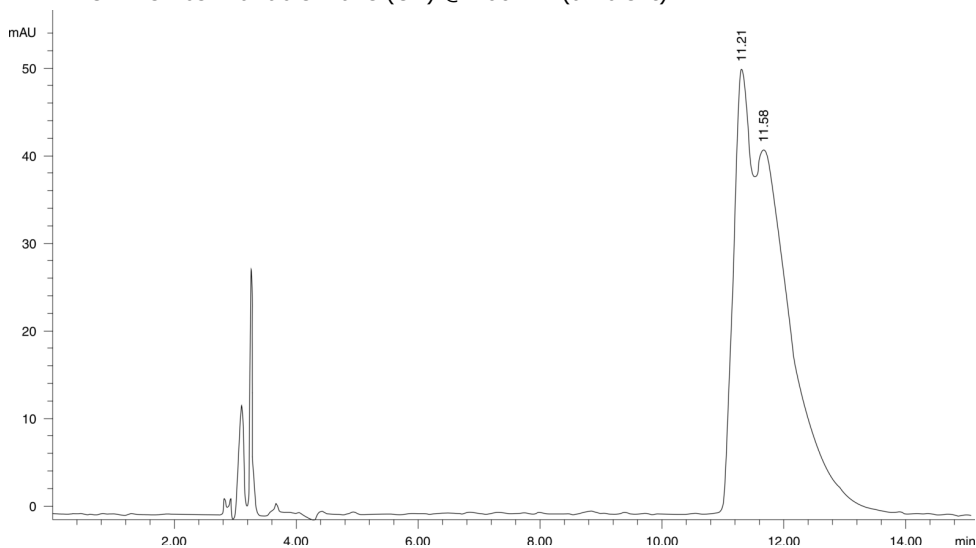