


Tetrahydrozoline

Column: Chirex® (S)-VAL and (R)-NEA, LC Column 250 x 4.6 mm, Ea
Dimensions: 250 x 4.6 mm ID
Order No: 00G-3014-E0
Elution Type: Isocratic
Eluent A: Hexane/Dichloroethane/Ethanol-TFA (55:35:10) w/ethanol-TFA (20:1)
Gradient


Step No.	Time (min)	Pct A
1	0	100

Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: UV-Vis Abs.-Variable Wave.(UV) @ 242 nm (ambient)

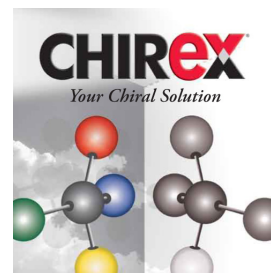

Products used in this application:

CHIREX®

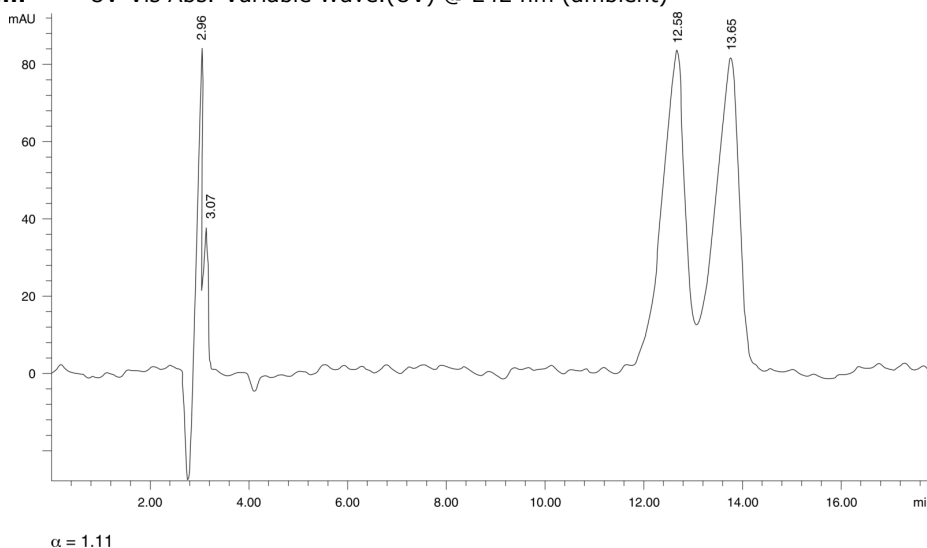

ANALYTES:

- 1 Tetrahydrozoline
Retention Time: 12.58 min
- 2 Tetrahydrozoline
Retention Time: 13.65 min

