


Thioridazine**Column:** N/A**Dimensions:** N/A**Order No:** N/A**Elution Type:** N/A**Eluent A:** Hexane/Dichloroethane/Ethanol-TFA (55:35:7) w/ ethanol-TFA (20:1)**Gradient** **Step No.** **Time (min)** **Pct A****Profile:** **1** 0 100**Flow Rate:** 1 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 264 nm (ambient)**ANALYTES:****1** Thioridazine

Retention Time: 11.7 min

2 Thioridazine

Retention Time: 12.61 min

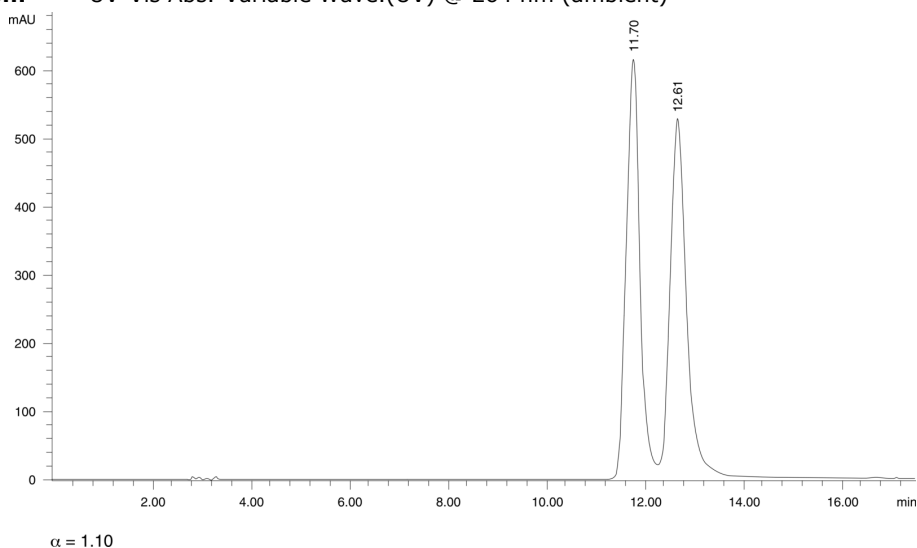