

Thioridazine

Column: Chirex® (S)-LEU and (R)-NEA, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-3020-E0

Elution Type: Isocratic


Eluent A: Hexane/Dichloroethane/Ethanol-TFA (77:20:3)

Gradient	Step No.	Time (min)	Pct A
Profile:	1	0	100

Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 264 nm (ambient)

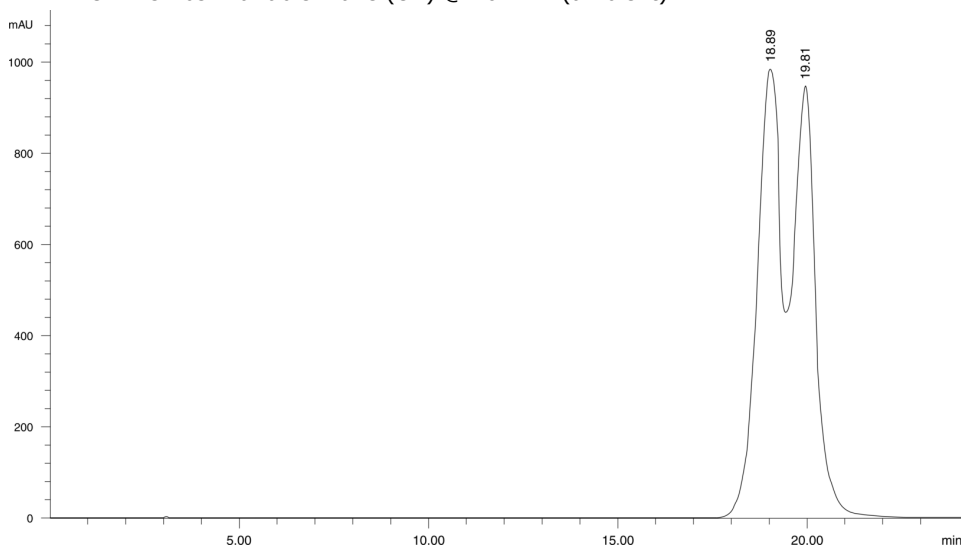

Products used in this application:

CHIREX®

ANALYTES:

- 1 Thioridazine
Retention Time: 18.89 min
- 2 Thioridazine
Retention Time: 19.81 min

