

Warfarin

Column: N/A

Dimensions: N/A

Order No: N/A

Elution Type: N/A

Eluent A: Hexane/Dichloroethane/Ethanol-TFA (88:10:2) ethanol-TFA 20:1


Gradient **Step No.** **Time (min)** **Pct A**

Profile: **1** 0 100

Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 282 nm (ambient)

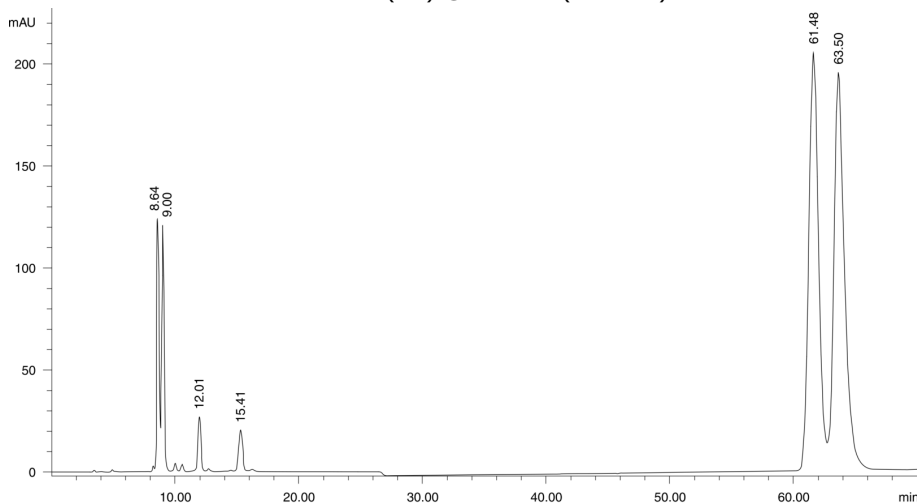

$\alpha = 1.03$

ANALYTES:

1 Warfarin

Retention Time: 61.48 min

2 Warfarin

Retention Time: 63.5 min

