


Tyrosine

Column: Chirex® 3126 (D)-penicillamine, LC Column 150 x 4.6 mm, Ea
Dimensions: 150 x 4.6 mm ID
Order No: 00F-3126-E0
Elution Type: Isocratic
Eluent A: 2 mM Copper (II) sulfate in water
Eluent B: Methanol
Gradient

Step No.	Time (min)	Pct A	Pct B
1	0	85	15


Profile:
Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

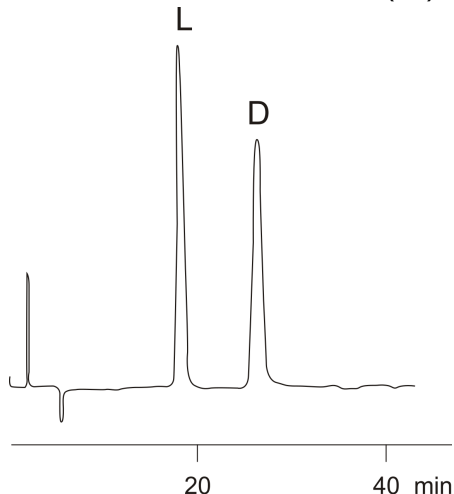

$\alpha = 1.34$

ANALYTES:

- 1 L-Tyrosine
- 2 D-Tyrosine

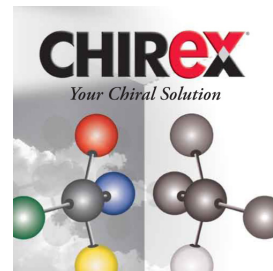

Products used in this application:

CHIREX

