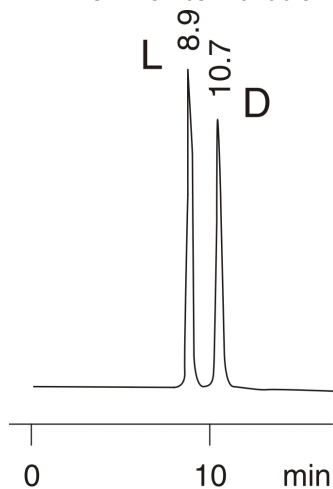


**N-3,5-Dinitrobenzoylvaline****Column:** N/A**Dimensions:** N/A**Order No:** N/A**Elution Type:** N/A**Eluent A:** Hexane**Eluent B:** Isopropanol**Eluent C:** Methanol**Eluent D:** Trifluoroacetic acid

| Gradient | Step No. | Time (min) | Pct A | Pct B | Pct C | Pct D |
|----------|----------|------------|-------|-------|-------|-------|
| Profile: | 1        | 0          | 70    | 20    | 10    | 0.2   |

**Flow Rate:** 1 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

$$\alpha = 1.30$$

**ANALYTES:**

- 1 L-N-3,5-Dinitrobenzoylvaline

Retention Time: 8.98 min

- 2 D-N-3,5-Dinitrobenzoylvaline

Retention Time: 10.65 min

