

p-Hydroxynorephedrine

Column: Chirex® 3126 (D)-penicillamine, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-3126-E0

Elution Type: Isocratic

Eluent A: p-Hydroxynorephedrine

Gradient


Step No.	Time (min)	Pct A
----------	------------	-------

Profile:	1	0	100
-----------------	----------	----------	------------


Flow Rate: 0.2 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 280 nm (ambient)

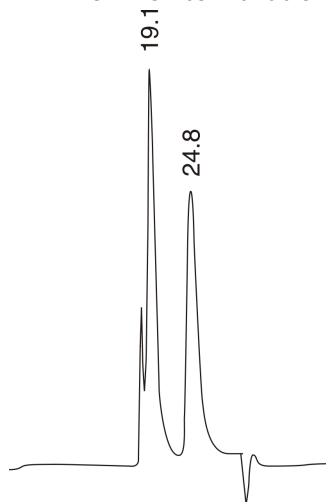

$$\alpha = 1.43$$

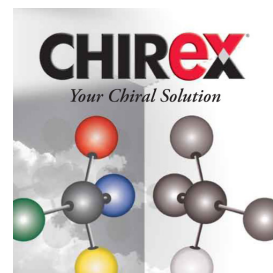

Products used in this application:

CHIREX®

ANALYTES:

1 p-Hydroxynorephedrine

Retention Time: 19.05 min

2 p-Hydroxynorephedrine

Retention Time: 24.76 min

