

2-Amino-n-butyric acid

Column: Chirex® 3126 (D)-penicillamine, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-3126-E0

Elution Type: Isocratic

Eluent A: 1 mM Copper (II) sulfate in water


Gradient	Step No.	Time (min)	Pct A
----------	----------	------------	-------

Profile:	1	0	100
----------	---	---	-----


Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

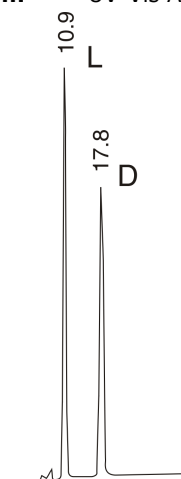

$$\alpha = 1.80$$

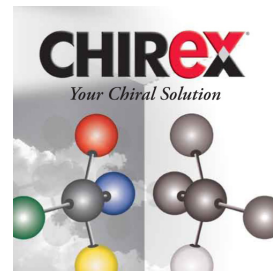

Products used in this application:

CHIREX®

ANALYTES:

- 1 L-2-Amino-n-butyric acid

Retention Time: 10.85 min

- 2 2-Amino-n-butyric acid

Retention Time: 17.703 min

