

Alloisoleucine

Column: Chirex® 3126 (D)-penicillamine, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-3126-E0

Elution Type: Isocratic

Eluent A: 2 mM Copper (II) sulfate in water


Eluent B: Methanol

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
1	0	85	15	

Flow Rate: 1 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

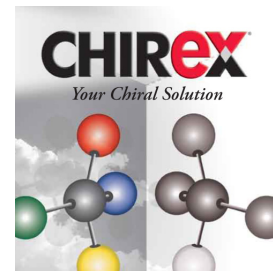

Products used in this application:

CHIREX

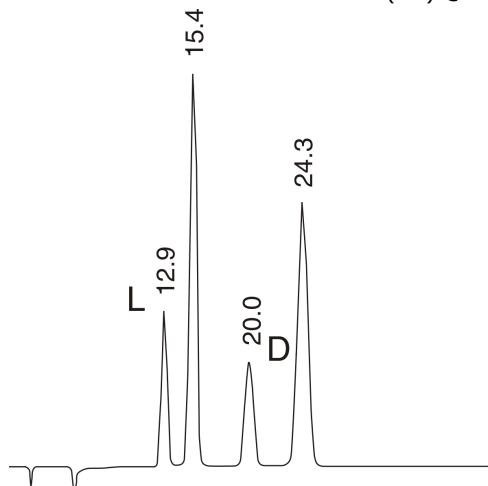

$$\alpha = 1.67$$

ANALYTES:

- 1** L-Alloisoleucine
Retention Time: 15.367 min
- 2** D-Alloisoleucine
Retention Time: 24.333 min

