

3-Amino-e-caprolactam

Column: Chirex® 3126 (D)-penicillamine, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-3126-E0

Elution Type: Isocratic



Eluent A: 1 mM Copper (II) sulfate in water

Gradient	Step No.	Time (min)	Pct A
Profile:	1	0	100

Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

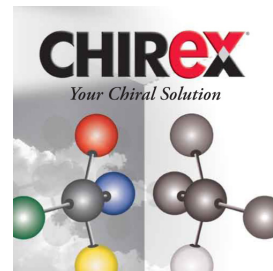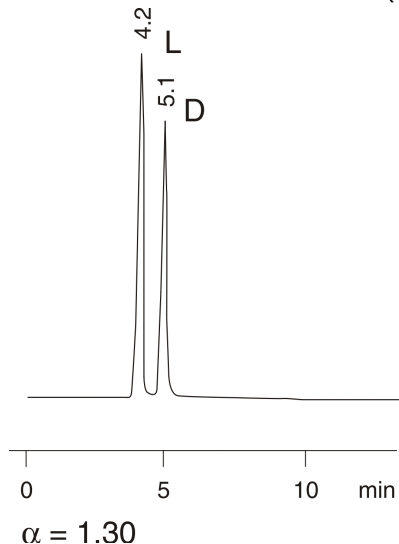

Products used in this application:

CHIREX

ANALYTES:

1 L-3-Amino-e-caprolactam

Retention Time: 4.21 min

2 Amino-e-caprolactam-d3

Retention Time: 5.05 min

