
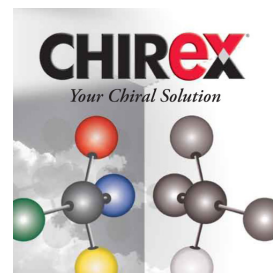


Leucylglycyl-glycine

Column: Chirex® 3126 (D)-penicillamine, LC Column 150 x 4.6 mm, Ea
Dimensions: 150 x 4.6 mm ID
Order No: 00F-3126-E0
Elution Type: Isocratic
Eluent A: 2 mM Copper (II) sulfate in water
Eluent B: Methanol
Gradient

Step No.	Time (min)	Pct A	Pct B
1	0	70	30

Profile:
Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

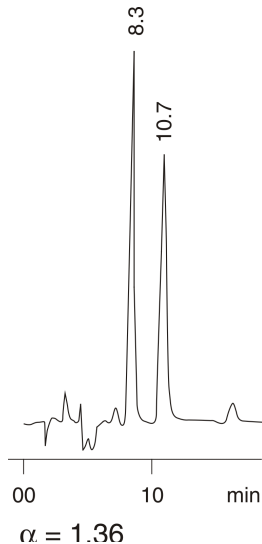

Products used in this application:

CHIREX

ANALYTES:

- 1 Leucylglycyl-glycine
Retention Time: 8.29 min
- 2 Leucylglycyl-glycine
Retention Time: 10.74 min

