

Tetrahydro-3-furancarboxylic acid

Column: Chirex® 3126 (D)-penicillamine, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-3126-E0

Elution Type: Isocratic


Eluent A: 1 mM Copper (II) sulfate in water

Gradient	Step No.	Time (min)	Pct A
Profile:	1	0	100


Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

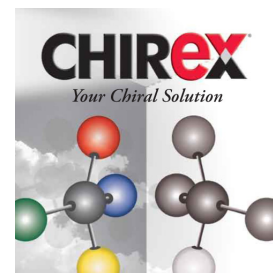

Products used in this application:

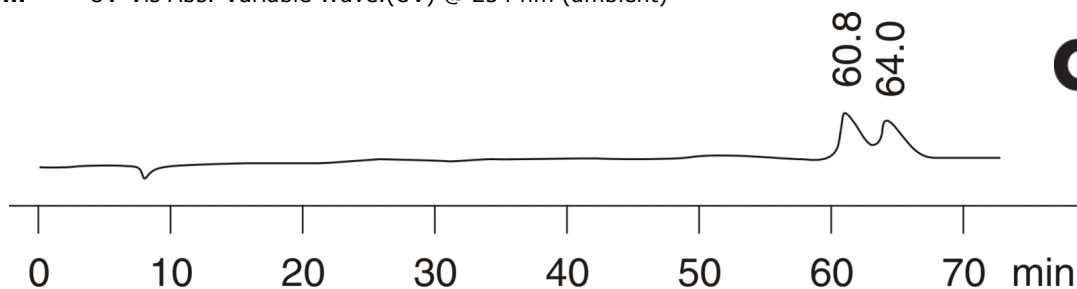The Chirex logo is displayed in a large, bold font.

$\alpha = 1.05$

ANALYTES:

- 1 Tetrahydro-3-furancarboxylic acid

Retention Time: 60.892 min

- 2 Tetrahydro-3-furancarboxylic acid

Retention Time: 64.021 min

