

Normetanephine

Column: Chirex® 3126 (D)-penicillamine, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-3126-E0

Elution Type: Isocratic

Eluent A: 5 mM Copper (II) acetate in water [pH 6.0]



Gradient

Step No.	Time (min)	Pct A
1	0	100

Flow Rate: 0.5 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 280 nm (ambient)

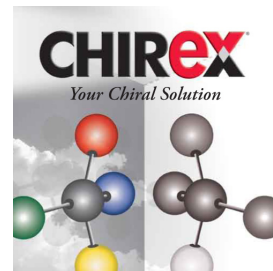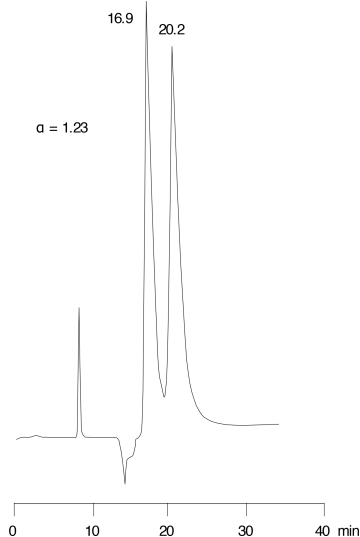

Products used in this application:

CHIREX®

ANALYTES:

1 Normetanephine

Retention Time: 16.9 min

2 Normetanephine

Retention Time: 20.2 min

