


Simultaneous Analysis of Organic Acids**Column:** N/A**Dimensions:** N/A**Order No:** N/A**Elution Type:** N/A**Eluent A:** 3.2 mM Na₂CO₃ / 1.0 mM NaHCO₃**Gradient** **Step No.** **Time (min)** **Pct A****Profile:** **1** 0 100**Flow Rate:** 0.7 mL/min**Col. Temp.:** 25 °C**Detection:** Conductivity @ (ambient)**Detector Info:** Metrohm 761 Compact IC

14176

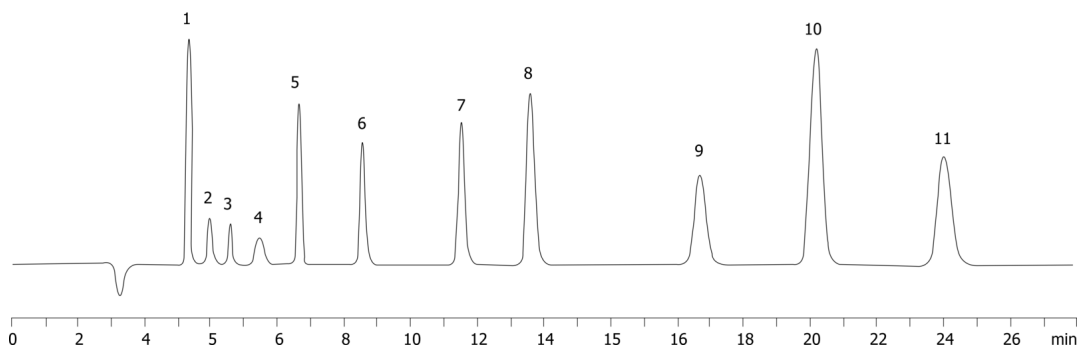

Simultaneous Analysis of Organic Acids

ANALYTES:

- 1 Fluoride
- 2 Acetate
- 3 Formate
- 4 Methacrylate
- 5 Chloride
- 6 Nitrite
- 7 Bromide
- 8 Nitrate
- 9 Phosphate
- 10 Sulfate
- 11 Oxalate

