

Beta-blockers on Synergi 4um Max-RP

Column: Synergi™ 4 µm Max-RP 80 Å, LC Column 50 x 4.6 mm, Ea

Dimensions: 50 x 4.6 mm ID

Order No: 00B-4337-E0

Elution Type: Gradient

Eluent A: 20mM KH₂PO₄, pH 2.5

Eluent B: Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	5	50	50

Flow Rate: 1.5 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 230 nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C12 4 x 3.0mm, 10/Pk Part No.: AJ0-6074

- Holder Part No.: KJ0-4282

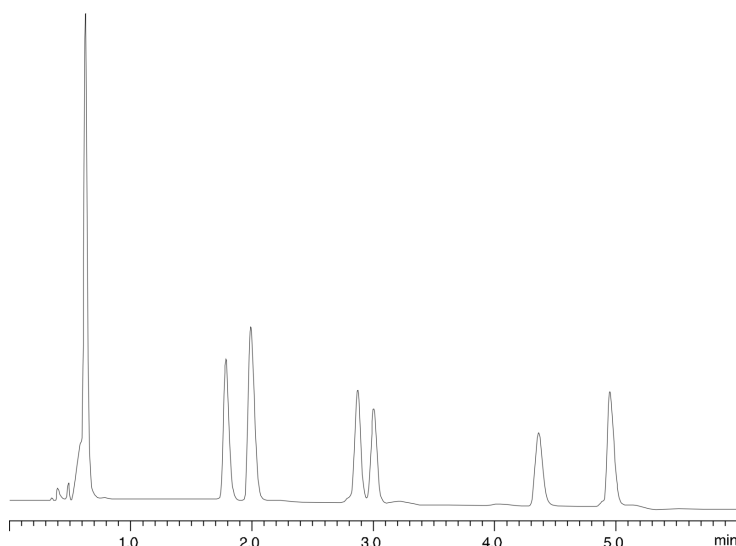

ANALYTES:

- 1 Atenolol
- 2 Nadolol
- 3 Pindolol
- 4 Acebutolol
- 5 Metoprolol
- 6 Labetalol
- 7 Alprenolol

