

## Lactoglobulin tryptic map on Luna C18(2) - 0.05% TFA

**Column:** Luna® 5 µm C18(2) 100 Å, LC Column 250 x 4.6 mm, Ea

**Dimensions:** 250 x 4.6 mm ID

**Order No:** 00G-4252-E0

**Elution Type:** Gradient

**Eluent A:** Water with 0.05% TFA

**Eluent B:** Acetonitrile with 0.05% TFA

| Gradient Profile: | Step No. | Time (min) | Pct A | Pct B |
|-------------------|----------|------------|-------|-------|
|                   | 1        | 0          | 97    | 3     |
|                   | 2        | 3          | 97    | 3     |
|                   | 3        | 90         | 60    | 40    |

**Flow Rate:** 1 mL/min


**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 210 nm (ambient)

**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C18 4 x 3.0mm, 10/Pk Part No.: AJ0-4287

- Holder Part No.: KJ0-4282

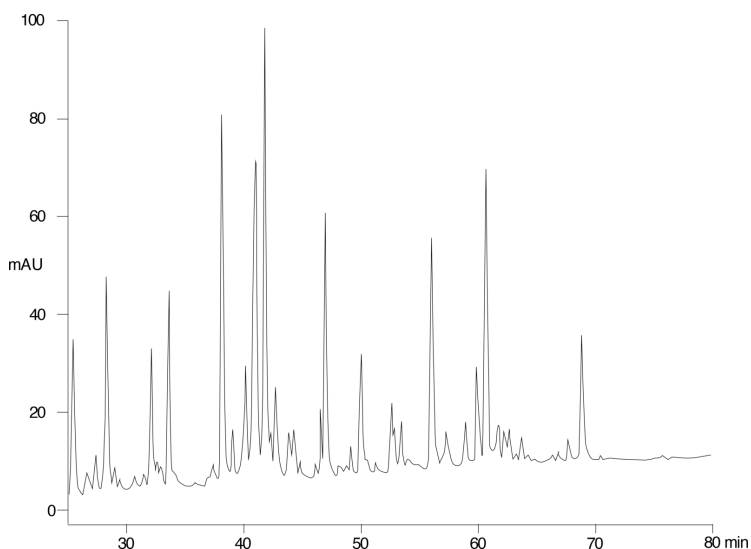

### ANALYTES:

1 Protein



Products used in this application:

