

Catecholamines on Synergi Hydro-RP

Column: Synergi™ 4 µm Hydro-RP 80 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4375-E0

Elution Type: Isocratic

Eluent A: 20mM Potassium phosphate pH 2.5

Gradient	Step No.	Time (min)	Pct A
Profile:	1	0	100

Flow Rate: 1 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 210 nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, AQ C18 4 x 3.0mm, 10/Pk Part No.: AJ0-7511

- Holder Part No.: KJ0-4282

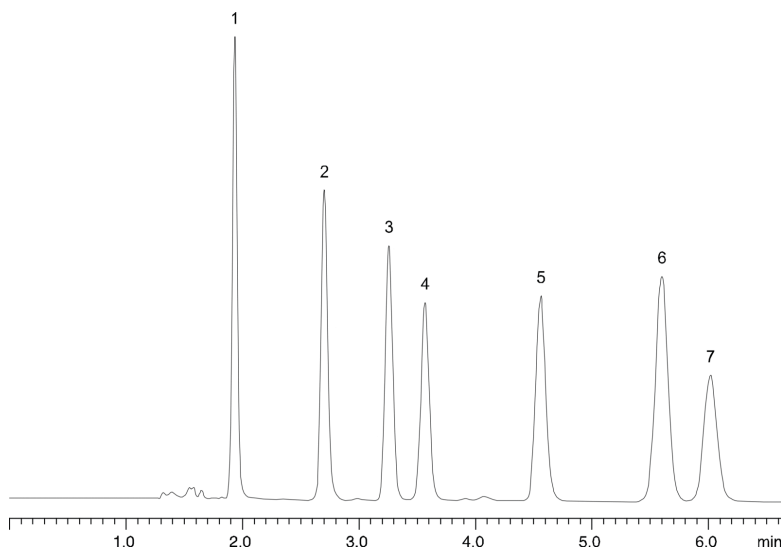

ANALYTES:

- 1 Norepinephrine
- 2 Epinephrine
- 3 6-Hydroxydopamine
- 4 Normetanephrine
- 5 Dopamine
- 6 L-DOPA
- 7 Epinine

