

Pesticides on Synergi Max-RP

Column: Synergi[™] 4 μ m Max-RP 80 Å, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4337-E0

Elution Type: Gradient

Eluent A: Water with 0.1% H₃PO₄

Eluent B: Acetonitrile with 0.1% H₃PO₄

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	30	5	95

Flow Rate: 1.5 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 210 nm (ambient)

Analyst Note: SecurityGuard[™] Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C12 4 x 3.0mm, 10/Pk Part No.: AJ0-6074

- Holder Part No.: KJ0-4282

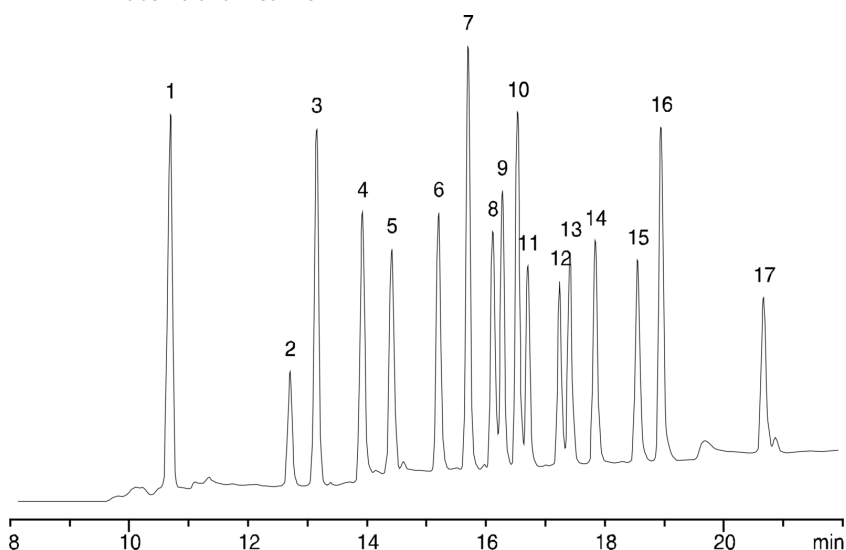

Pesticides on Synergi Max-RP

ANALYTES:

- 1** Desethylatrazine
- 2** Metoxuron
- 3** Hexazinone
- 4** Simazine
- 5** Cyanazine
- 6** Methabenzthiazuron
- 7** Atrazine
- 8** Chlorotoluron
- 9** Monolinuron
- 10** Isoproturon
- 11** Diuron
- 12** Metobromuron
- 13** Metazachlor
- 14** Sebuthylazine
- 15** Terbutylazine
- 16** Linuron
- 17** Metolachlor

