

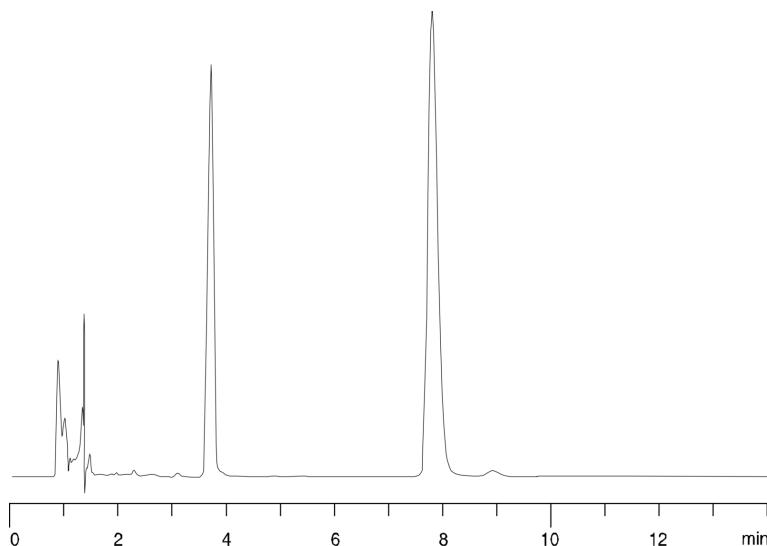
Goldenseal on Synergi Hydro-RP**Column:** SynergiTM 4 μ m Hydro-RP 80 Å, LC Column 150 x 4.6 mm, Ea**Dimensions:** 150 x 4.6 mm ID**Order No:** 00F-4375-E0**Elution Type:** Isocratic**Eluent A:** 100mM Dihydrogen Potassium phosphate, pH 4.6**Eluent B:** Acetonitrile

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	73	27

Flow Rate: 1.5 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 235 nm (ambient)**Analyst Note:** SecurityGuardTM Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, AQ C18 4 x 3.0mm, 10/Pk Part No.: AJ0-7511

- Holder Part No.: KJ0-4282

**ANALYTES:**

- 1 Hydrastine
- 2 Berberine

