

Analgesics on Synergi Hydro-RP

Column: Synergi™ 4 µm Hydro-RP 80 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4375-E0

Elution Type: Isocratic

Eluent A: Acetonitrile

Eluent B: 20mM Potassium phosphate, pH 2.5

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
1	0	50	50	

Flow Rate: 1.5 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 230 nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, AQ C18 4 x 3.0mm, 10/Pk Part No.: AJ0-7511

- Holder Part No.: KJ0-4282

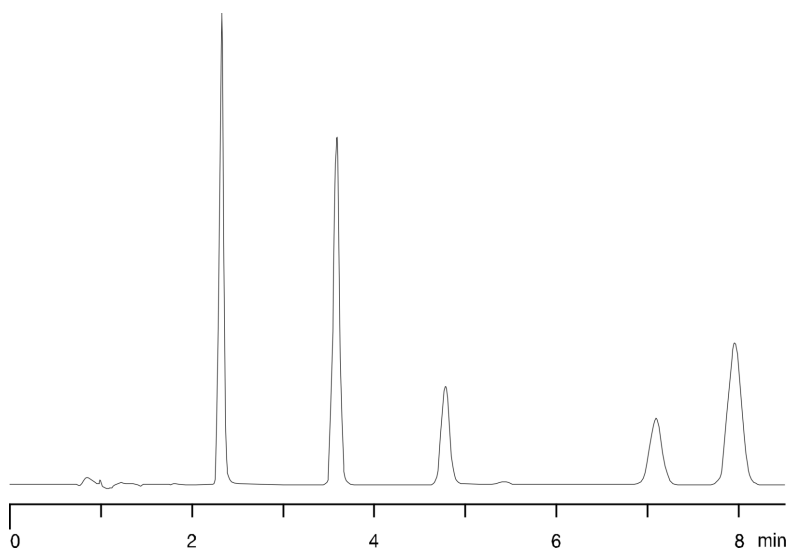

ANALYTES:

- 1 Indoprofen / Ethyl paraben
- 2 Naproxen
- 3 Diflunisal
- 4 Indomethacin
- 5 Ibuprofen

