

Antidepressants on Luna C8(2)

Column: Luna® 5 µm C8(2) 100 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4249-E0

Elution Type: Isocratic

Eluent A: 20mM Potassium phosphate, pH 7

Eluent B: Acetonitrile

Eluent C: Methanol

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B	Pct C
1	1	0	30	35	35

Flow Rate: 1.5 mL/min


Col. Temp.: 30 °C

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C8 4 x 3.0mm, 10/Pk Part No.: AJ0-4290

- Holder Part No.: KJ0-4282

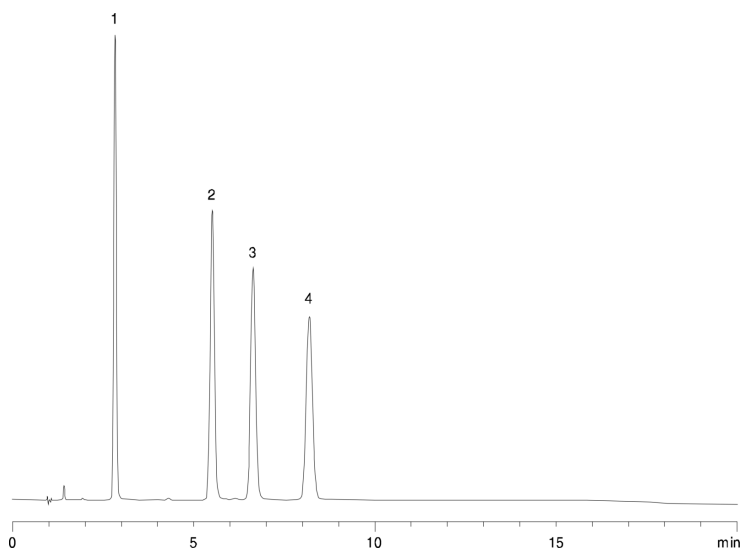

ANALYTES:

- 1 Desipramine
- 2 Nortriptyline
- 3 Amitriptyline
- 4 Clomipramine



Products used in this application:

