

Simple Sugars (#4) on Luna Amino

Column: Luna® 5 µm NH2 100 Å, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4378-E0

Elution Type: Isocratic

Eluent A: Acetonitrile

Eluent B: Water

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	80	20

Flow Rate: 3 mL/min

Col. Temp.: 40 °C

Detection: Refractive Index (RI) @ (ambient)

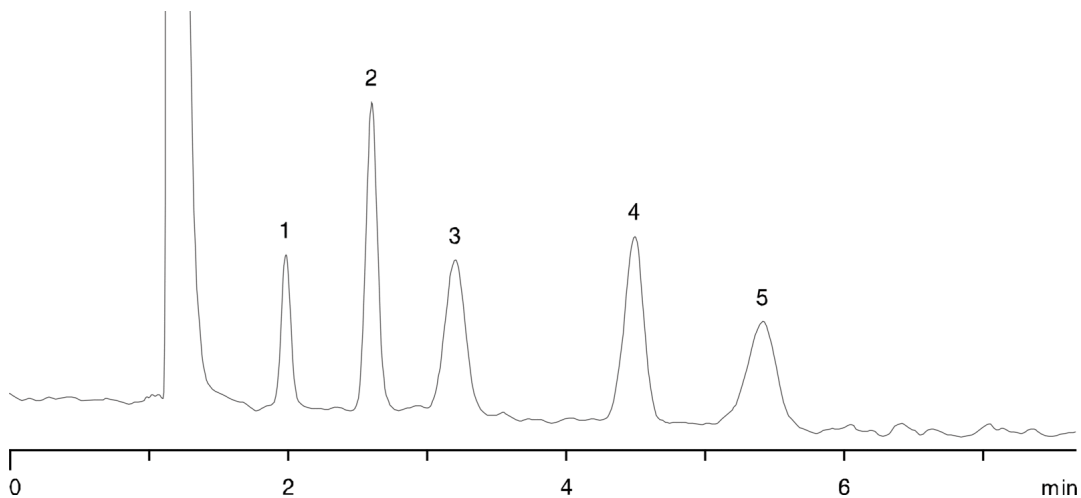
Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, NH2 4 x 3.0mm, 10/Pk Part No.: AJ0-4302

- Holder Part No.: KJ0-4282



Products used in this application:

**ANALYTES:**

- 1 Rhamnose
- 2 Fructose
- 3 Galactose
- 4 Sucrose
- 5 Maltose