

K vitamins on Luna Amino

Column: Luna® 5 µm NH2 100 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4378-E0

Elution Type: Isocratic

Eluent A: Hexane

Gradient Profile:	Step No.	Time (min)	Pct A
	1	0	100

Flow Rate: 1.5 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 270 nm (ambient)


Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, NH2 4 x 3.0mm, 10/Pk Part No.: AJ0-4302

- Holder Part No.: KJ0-4282



Products used in this application:

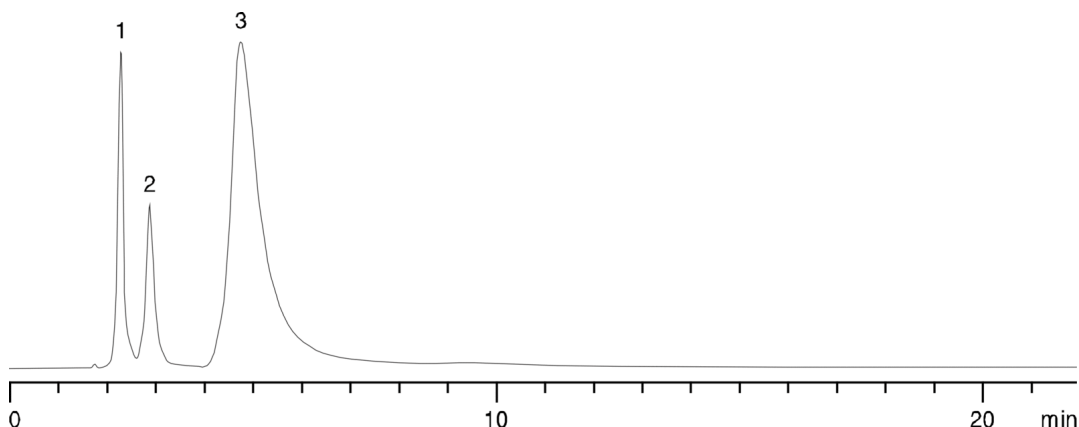

ANALYTES:

- 1 Vitamin K1
- 2 Vitamin K2
- 3 Vitamin K3

