


EPA Method 624 CLP

Column: Zebron™ ZB-624, GC Cap. Column 30 m x 0.53 mm x 3.00 µm, Ea**Phase:** 6% Cyanopropylphenyl 94% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.53 mm x 3 µm**Order No:** 7HK-G005-36**Oven Profile:** 35°C (6 min), 5°C/min to 140°C (8 min)**Carrier Gas:** Constant Flow Helium, 2.65 mL/min**Injection:** Split 10:1 1 µL @ 250°C**Detection:** Flame Ionization (FID) (350°C)

Products used in this application:

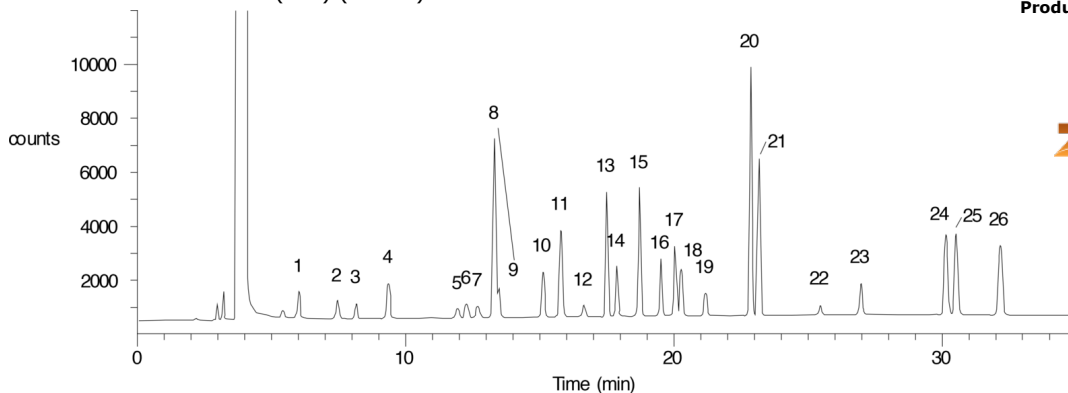

EPA Method 624 CLP

ANALYTES:

- 1** 1,1-Dichloroethene
- 2** Methylene chloride
- 3** trans-1,2-Dichloroethene
- 4** 1,1-Dichloroethane
- 5** Chloroform
- 6** 1,1,1-Trichloroethane
- 7** Carbon tetrachloride
- 8** Benzene
- 9** 1,2-Dichloroethane
- 10** Trichloroethene
- 11** 1,2-Dichloropropane
- 12** Bromodichloromethane
- 13** 2-Chloroethyl vinyl ether
- 14** cis-1,3-Dichloropropene
- 15** Toluene
- 16** trans-1,3-Dichloropropene
- 17** 1,1,2-Trichloroethane
- 18** Tetrachloroethane
- 19** Dibromochloromethane
- 20** Ethylbenzene
- 21** Chlorobenzene
- 22** Bromoform
- 23** 1,1,2,2-Tetrachloroethane
- 24** 1,3-Dichlorobenzene
- 25** 1,4-Dichlorobenzene
- 26** 1,2-Dichlorobenzene

