

## TFA Concentration Comparison (6) - Tryptic Digest of Beta-Lactoglobulin

**Column:** Jupiter® 4 µm Proteo 90 Å, LC Column 250 x 4.6 mm, Ea

**Dimensions:** 250 x 4.6 mm ID

**Order No:** 00G-4396-E0

**Elution Type:** Gradient

**Eluent A:** 1.25mM TFA

**Eluent B:** 1.25mM TFA in Acetonitrile

| Gradient Profile: | Step No. | Time (min) | Pct A | Pct B |
|-------------------|----------|------------|-------|-------|
|                   | 1        | 0          | 95    | 5     |
|                   | 2        | 5          | 95    | 5     |
|                   | 3        | 60         | 55    | 45    |

**Flow Rate:** 1 mL/min


**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 210 nm (ambient)


**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C12 4 x 3.0mm, 10/Pk Part No.: AJ0-6074

- Holder Part No.: KJ0-4282

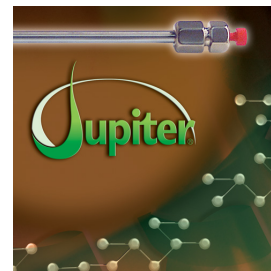

Products used in this application:

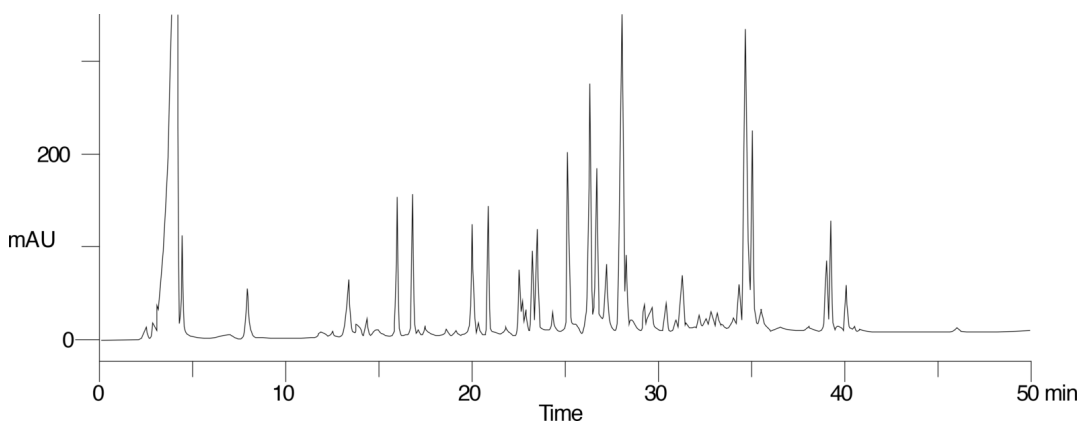

### ANALYTES:

1 beta-Lactoglobulin

