

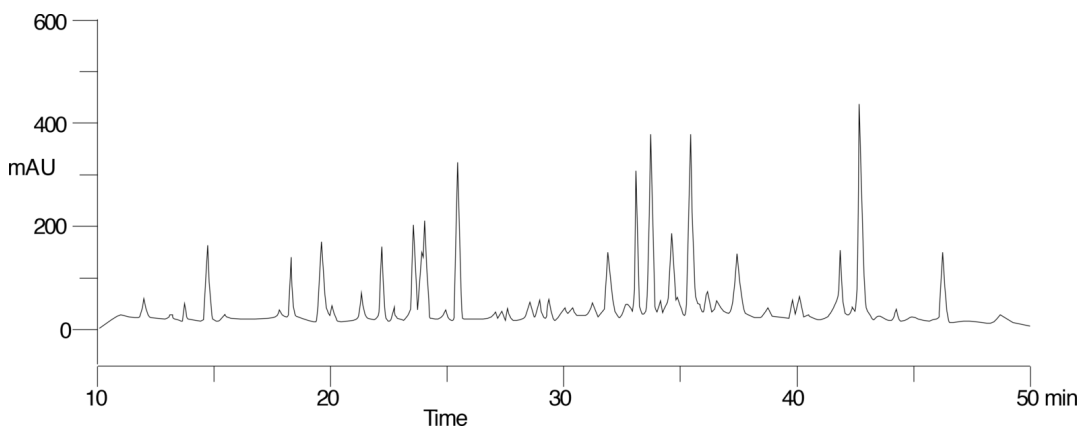
Tryptic Digest of beta-Lactoglobulin (oxidation) on Jupiter Proteo**Column:** Jupiter® 4 µm Proteo 90 Å, LC Column 250 x 4.6 mm, Ea**Dimensions:** 250 x 4.6 mm ID**Order No:** 00G-4396-E0**Elution Type:** Gradient**Eluent A:** 1.25mM TFA in Water**Eluent B:** 1.0mM TFA in Acetonitrile

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	95	5
	2	5	95	5
	3	60	60	40

Flow Rate: 1 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 210 nm (ambient)**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C12 4 x 3.0mm, 10/Pk Part No.: AJ0-6074

- Holder Part No.: KJ0-4282

**Products used in this application:****ANALYTES:****1** beta-Lactoglobulin