

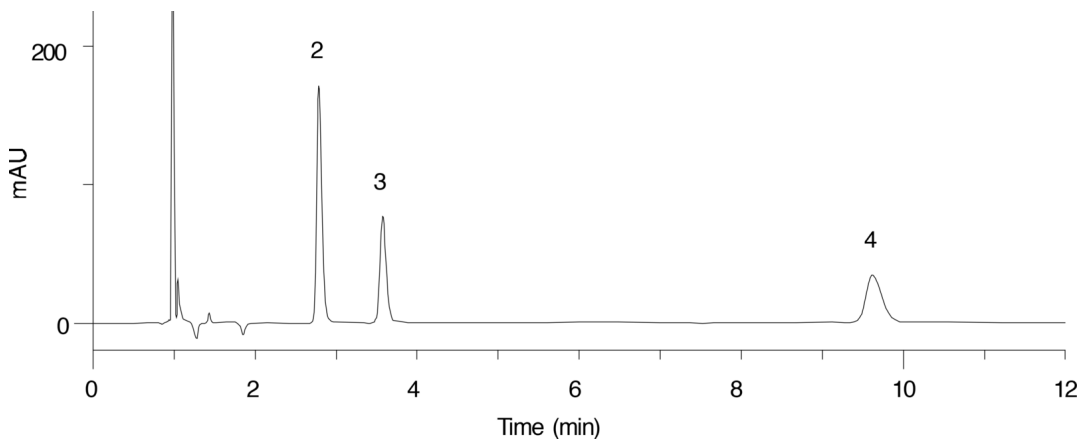
Cough and Colds on Synergi Hydro-RP**Column:** Synergi™ 4 µm Hydro-RP 80 Å, LC Column 150 x 4.6 mm, Ea**Dimensions:** 150 x 4.6 mm ID**Order No:** 00F-4375-E0**Elution Type:** Isocratic**Eluent A:** Potassium Phosphate Buffer pH 2.5**Eluent B:** Acetonitrile

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	75	25

Flow Rate: 1.5 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 210 nm (ambient)**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, AQ C18 4 x 3.0mm, 10/Pk Part No.: AJ0-7511

- Holder Part No.: KJ0-4282

**ANALYTES:**

- 1 Maleic acid
- 2 Triprolidine
- 3 Chlorpheniramine
- 4 Diphenhydramine

