

Carotenoids from Haematococcus Algae Powder

Column: Luna® 3 µm Silica (2) 100 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4162-E0

Elution Type: Isocratic

Eluent A: Hexane/Acetone (82:18, v/v)

Gradient Profile:	Step No.	Time (min)	Pct A
1	0	100	

Flow Rate: 1.2 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Fixed Wave. (UV) @ 474 nm (ambient)


Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Silica 4 x 3.0mm, 10/Pk Part No.: AJ0-4348

- Holder Part No.: KJ0-4282



Products used in this application:

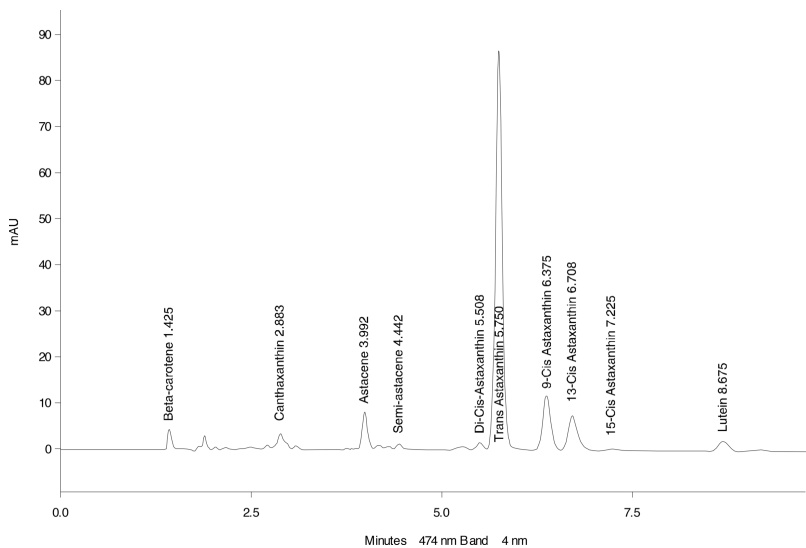

Carotenoids from Haematococcus Algae Powder

ANALYTES:

- 1** beta-Carotene
- 2** Canthaxanthin
- 3** Astacene
- 4** Semi-Astacene
- 5** Di-cis-astaxanthin
- 6** trans-Astaxanthin
- 7** 9-cis-Astaxanthin
- 8** 13-cis-Astaxanthin
- 9** 15-cis-Astaxanthin
- 10** Lutein

