

Extraction of Bronchodilators on Strata™ X

Column: Luna® 5 µm C8(2) 100 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4249-E0

Elution Type: Isocratic

Eluent A: 20 mM KH₂PO₄, pH 2.5

Eluent B: ACN

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5

Flow Rate: 1.5 mL/min


Col. Temp.: 25 °C

Detection: UV-Vis Abs.-Fixed Wave. (UV) @ 273 nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C8 4 x 3.0mm, 10/Pk Part No.: AJ0-4290

- Holder Part No.: KJ0-4282

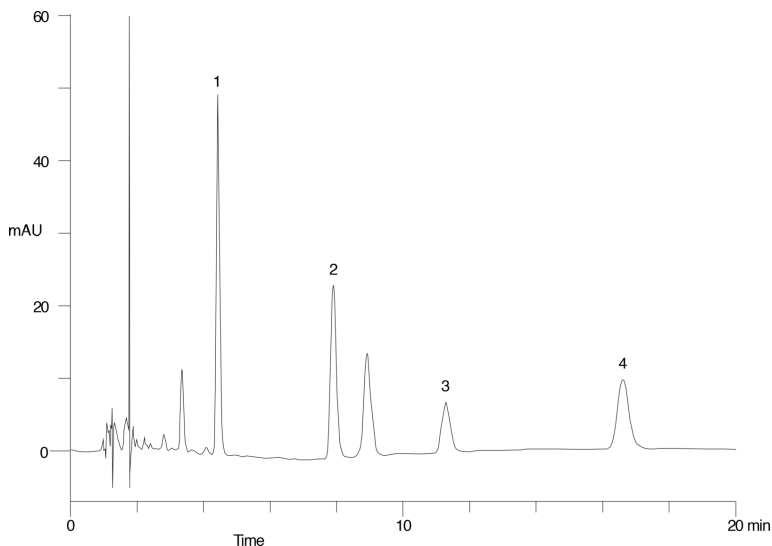

ANALYTES:

- 1 Theobromine
- 2 Theophylline
- 3 Beta hydroxy theophylline
- 4 Caffeine



Products used in this application:

