

## Deoxynucleosides on Luna SCX

**Column:** Luna® 5 µm SCX 100 Å, LC Column 150 x 4.6 mm, Ea

**Dimensions:** 150 x 4.6 mm ID

**Order No:** 00F-4398-E0

**Elution Type:** Isocratic

**Eluent A:** 50mM Potassium Phosphate, pH 2.5

| Gradient | Step No. | Time (min) | Pct A |
|----------|----------|------------|-------|
| Profile: | 1        | 0          | 100   |

**Flow Rate:** 1.5 mL/min

**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 210 nm (ambient)

**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, SCX 4 x 3.0mm, 10/Pk Part No.: AJ0-4308


- Holder Part No.: KJ0-4282



Products used in this application:



14557

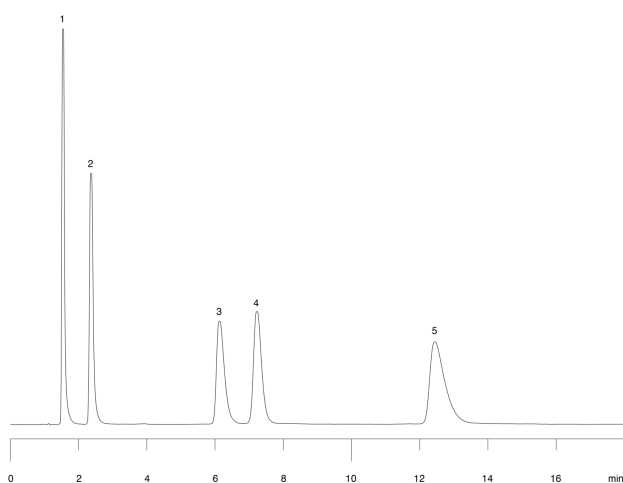

### ANALYTES:

- 1 d-Uridine
- 2 d-Inosine
- 3 d-Guanosine
- 4 d-Cytidine
- 5 d-Adenosine

