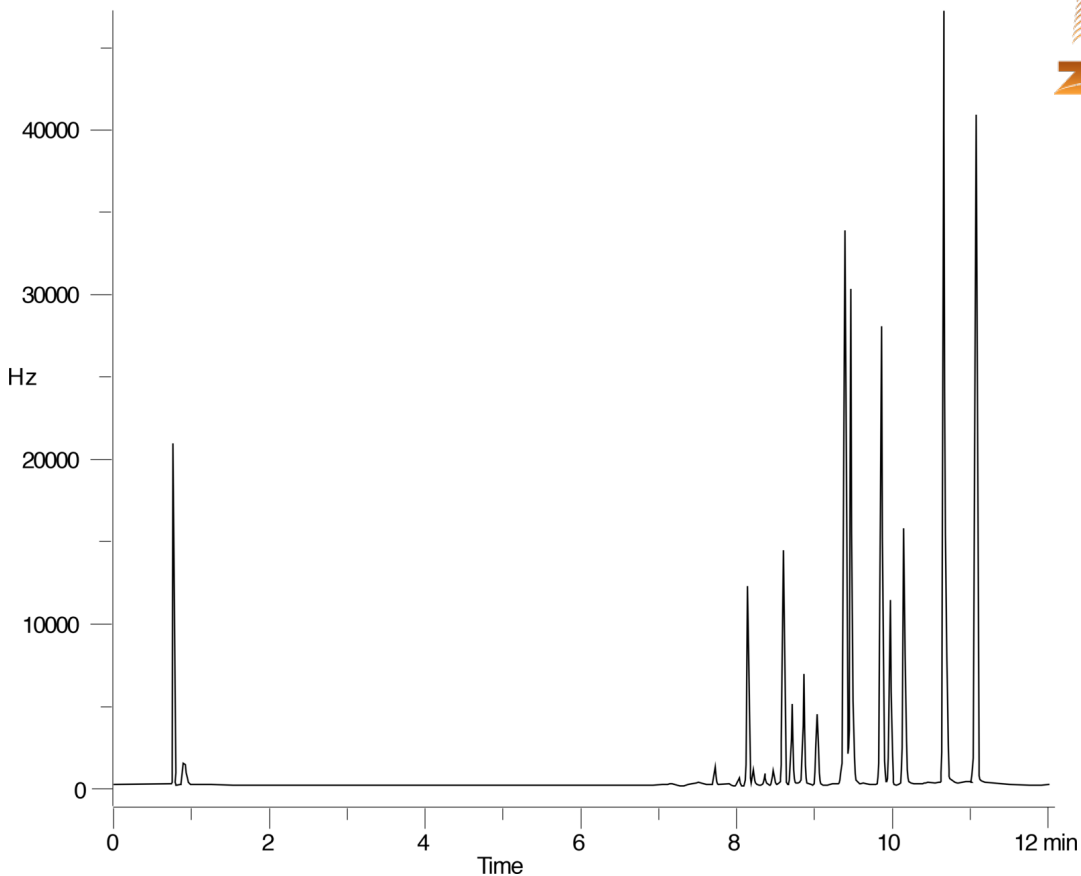


**Aroclor 1268 on ZB5 15m x .32mm x .25um****Column:** Zebtron<sup>™</sup> ZB-5, GC Cap. Column 15 m x 0.32 mm x 0.25  $\mu$ m, Ea**Phase:** 5% Phenyl 95% Dimethylpolysiloxane**Dimensions:** 15 meters x 0.32 mm x 0.25  $\mu$ m**Order No:** 7EM-G002-11**Oven Profile:** 150°C to 300°C at 12°C/min**Carrier Gas:** Constant Flow Helium, 1.7 mL/min**Injection:** split 15:1:1 1  $\mu$ L @ 250°C**Detection:** Electron Capture (ECD) (375°C)  
14576**Products used in this application:****ANALYTES:****1** Aroclor 1268