


Aroclors 1221 and 1268 on ZB5 15m x .32mm x .25um**Column:** Zebtron[™] ZB-5, GC Cap. Column 15 m x 0.32 mm x 0.25 μ m, Ea**Phase:** 5% Phenyl 95% Dimethylpolysiloxane**Dimensions:** 15 meters x 0.32 mm x 0.25 μ m**Order No:** 7EM-G002-11**Oven Profile:** 150°C to 300°C at 12°C/min**Carrier Gas:** Constant Flow Helium, 1.7 mL/min**Injection:** split 15:1:1 1 μ L @ 250°C**Detection:** Electron Capture (ECD) (375°C)
14579**Products used in this application:****ANALYTES:**

- 1 Aroclor 1221
- 2 Aroclor 1268

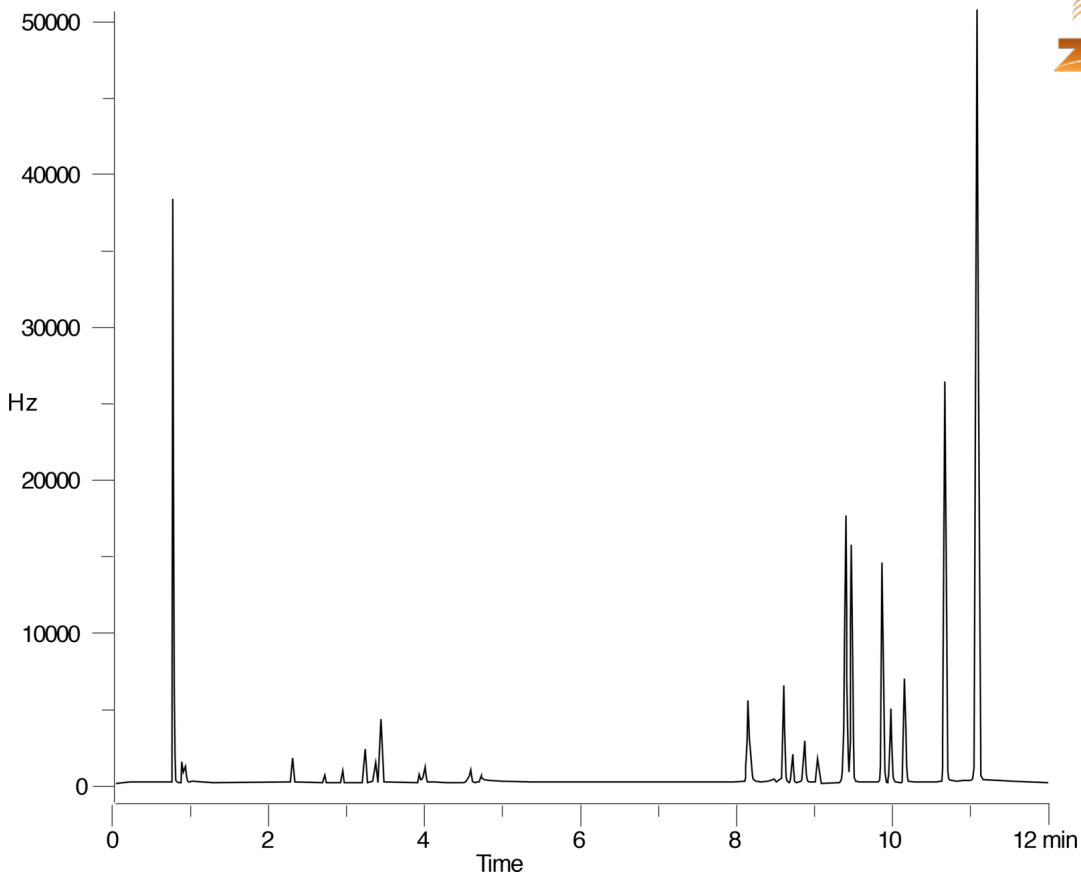