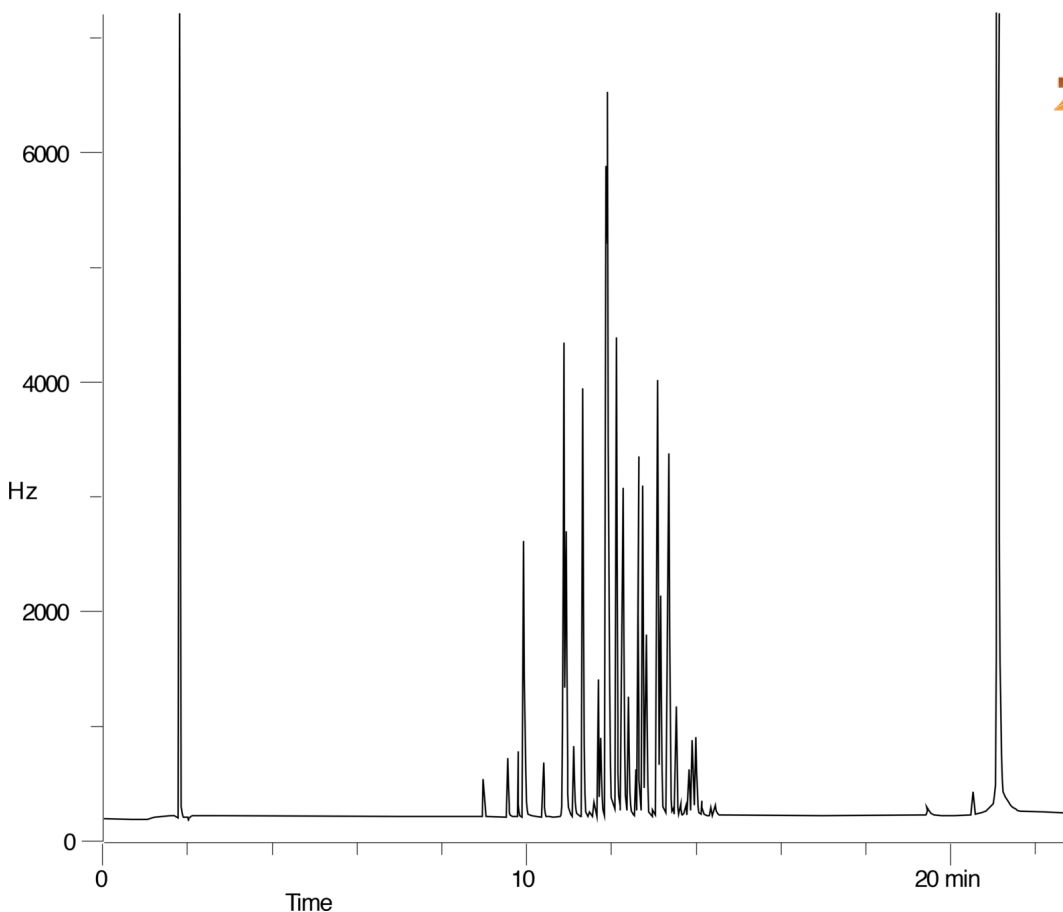


Aroclor 1016 on ZB5 30m x .32mm x .25um**Column:** Zebtron[™] ZB-5, GC Cap. Column 30 m x 0.32 mm x 0.25 μ m, Ea**Phase:** 5% Phenyl 95% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.32 mm x 0.25 μ m**Order No:** 7HM-G002-11**Oven Profile:** 120°C (1 min) to 300°C at 9°C/min (hold 5 min)**Carrier Gas:** Constant Flow Helium, 1.7 mL/min**Injection:** split 15:1:1 1 μ L @ 250°C**Detection:** Electron Capture (ECD) (375°C)
14580**Products used in this application:****ANALYTES:****1** Aroclor 1016