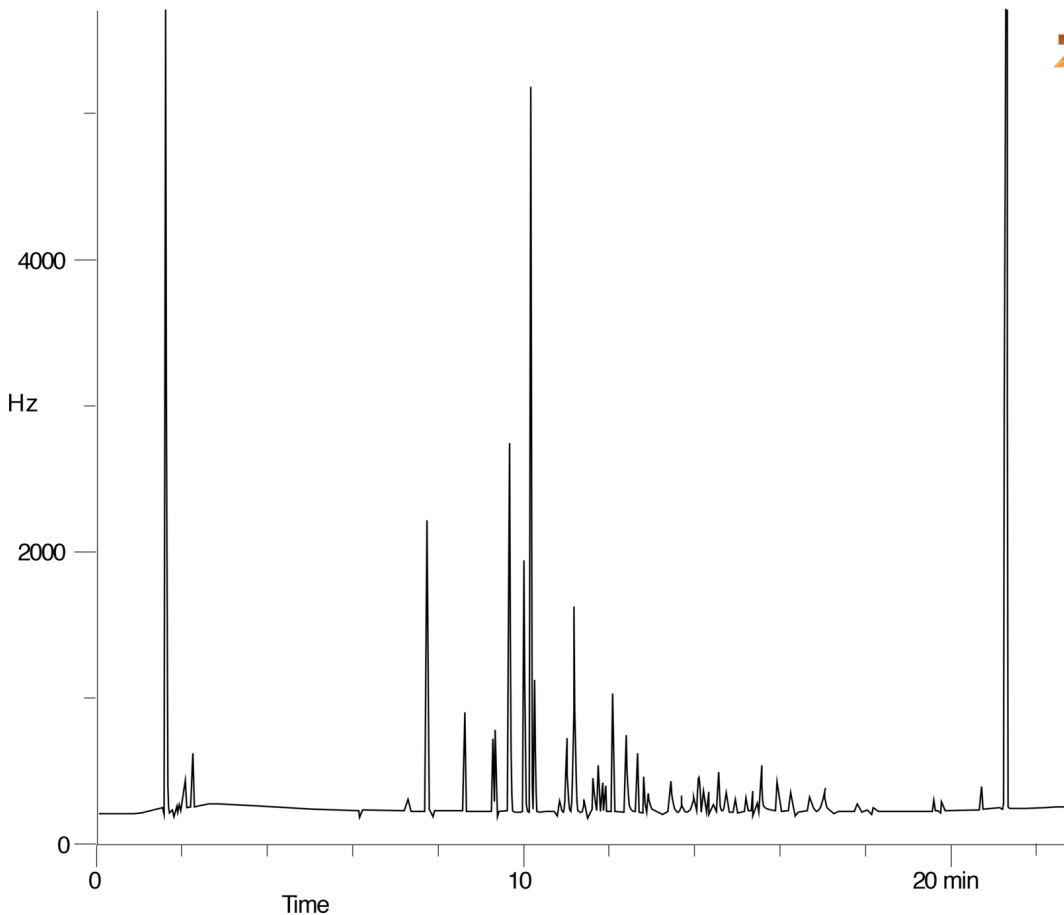


Aroclor 1221 on ZB35 30m x .32mm x .25um**Column:** Zebtron™ ZB-35, GC Cap. Column 30 m x 0.32 mm x 0.25 µm, Ea**Phase:** 35% Phenyl 65% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.32 mm x 0.25 µm**Order No:** 7HM-G003-11**Oven Profile:** 120°C (1 min) to 300°C at 9°C/min (hold 5 min)**Carrier Gas:** Constant Flow Helium, 1.7 mL/min**Injection:** Split :1 1 µL @ 250°C**Detection:** Electron Capture (ECD) (375°C)
14593**Products used in this application:****ANALYTES:****1** Aroclor 1221