


Aroclors 1016 + 1262 on ZB35 30m x .32mm x .25um**Column:** Zebtron[™] ZB-35, GC Cap. Column 30 m x 0.32 mm x 0.25 μ m, Ea**Phase:** 35% Phenyl 65% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.32 mm x 0.25 μ m**Order No:** 7HM-G003-11**Oven Profile:** 120°C (1 min) to 300°C at 9°C/min (hold 5 min)**Carrier Gas:** Constant Flow Helium, 1.7 mL/min**Injection:** Split 15:1:1 1 μ L @ 250°C**Detection:** Electron Capture (ECD) (375°C)
14602**Products used in this application:****ANALYTES:**

- 1 Aroclor 1016
- 2 Aroclor 1262

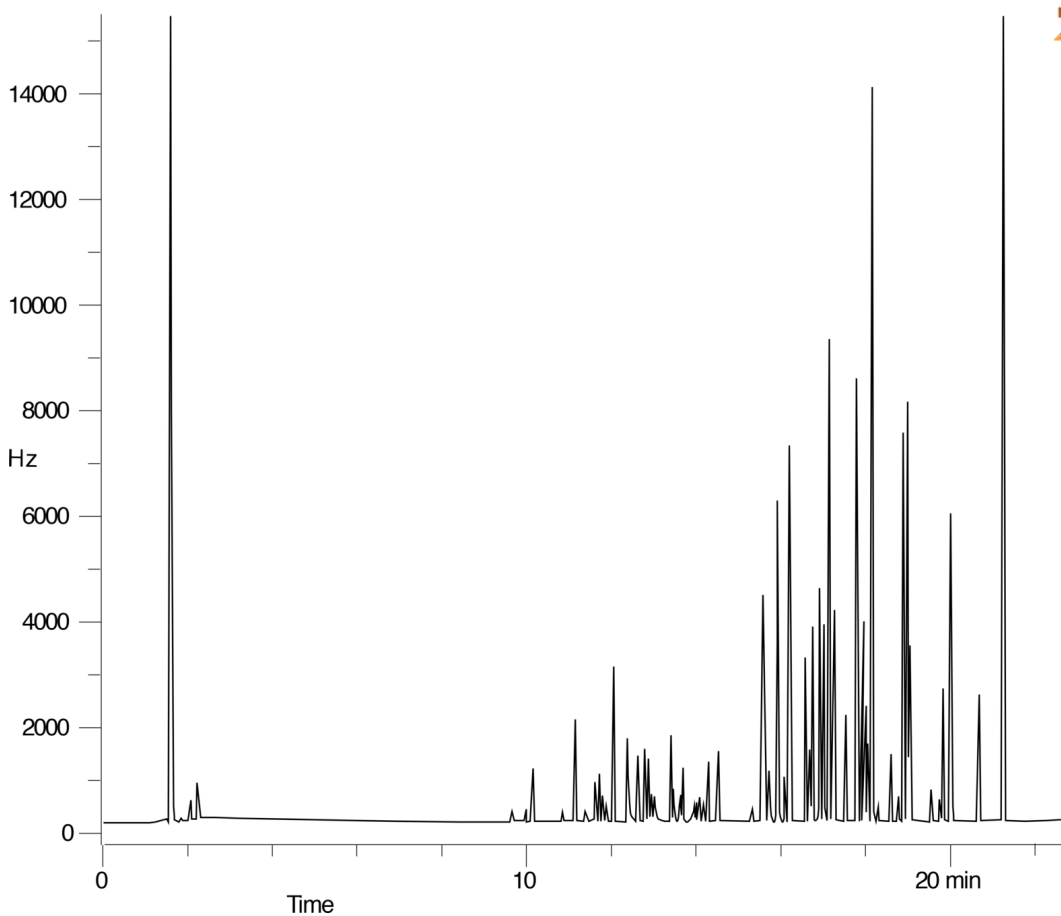