

Peptides from Serum on Strata™ X

Column: Jupiter® 4 µm Proteo 90 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4396-E0

Elution Type: Gradient

Eluent A: 0.1% TFA

Eluent B: ACN:0.1% TFA in water (90:10)

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	85	15
	2	22	40	60

Flow Rate: 1 mL/min


Col. Temp.: 30 °C

Detection: UV-Vis Abs.-Diode Array (PDA) @ 214 nm (ambient)


Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C12 4 x 3.0mm, 10/Pk Part No.: AJ0-6074

- Holder Part No.: KJ0-4282

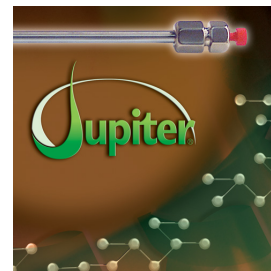

Products used in this application:

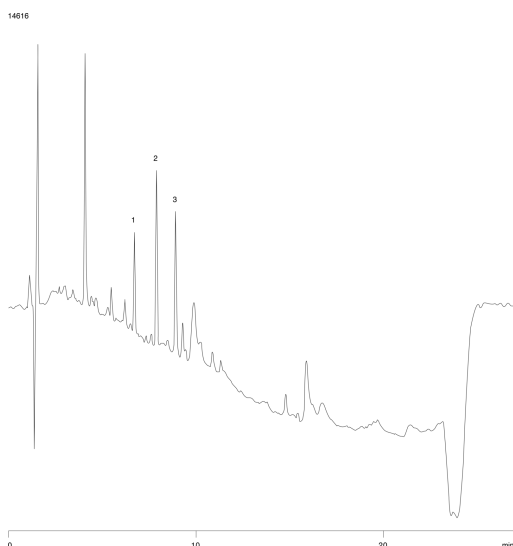

ANALYTES:

- 1** Oxytocin
- 2** Angiotensin II
- 3** Angiotensin I

