

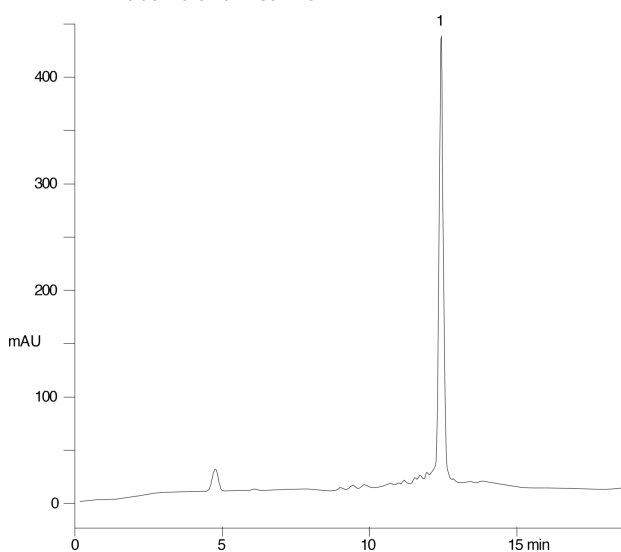
**20mer Hetero-Oligonucleotide on Jupiter Proteo****Column:** N/A**Dimensions:** N/A**Order No:** N/A**Elution Type:** N/A**Eluent A:** 0.1M TEAA, pH 7.0**Eluent B:** Acetonitrile

| Gradient        | Step No. | Time (min) | Pct A | Pct B |
|-----------------|----------|------------|-------|-------|
| <b>Profile:</b> | <b>1</b> | 0          | 97    | 3     |
|                 | <b>2</b> | 20         | 85    | 15    |

**Flow Rate:** 1 mL/min**Col. Temp.:** 60 °C**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 260 nm (ambient)**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C12 4 x 3.0mm, 10/Pk Part No.: AJ0-6074

- Holder Part No.: KJ0-4282

**ANALYTES:****1** 20mer heteroligonucleotide