

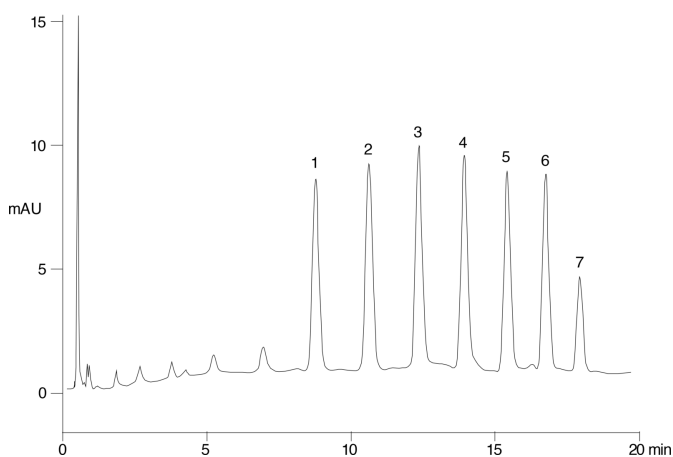
12-18mer Poly-dT Oligonucleotide Mix On Jupiter Proteo**Column:** N/A**Dimensions:** N/A**Order No:** N/A**Elution Type:** N/A**Eluent A:** 0.1M TEAA, pH 7.0**Eluent B:** Acetonitrile

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	97	3
	2	20	85	15

Flow Rate: 1 mL/min**Col. Temp.:** 60 °C**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 260 nm (ambient)**Analyst Note:** SecurityGuardTM Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C12 4 x 3.0mm, 10/Pk Part No.: AJ0-6074

- Holder Part No.: KJ0-4282

**ANALYTES:**

- 1 12-18mer Poly-dT oligonucleotide mixture

