

Anion Exchange 20mer Oligo on Luna Amino

Column: Luna® 5 µm NH₂ 100 Å, LC Column 150 x 4.6 mm, Ea
Dimensions: 150 x 4.6 mm ID
Order No: 00F-4378-E0
Elution Type: Gradient
Eluent A: 20mM Sodium Phosphate in 90/10 Water/Acetonitrile
Eluent B: 20mM Sodium Phosphate in 90/10 Water/Acetonitrile + 1000mM NaCl
Gradient

Step No.	Time (min)	Pct A	Pct B
1	2	75	25
2	45	5	95

Flow Rate: 1 mL/min


Col. Temp.: 40 °C

Detection: UV-Vis Abs.-Diode Array (PDA) @ nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, NH₂ 4 x 3.0mm, 10/Pk Part No.: AJ0-4302

- Holder Part No.: KJ0-4282

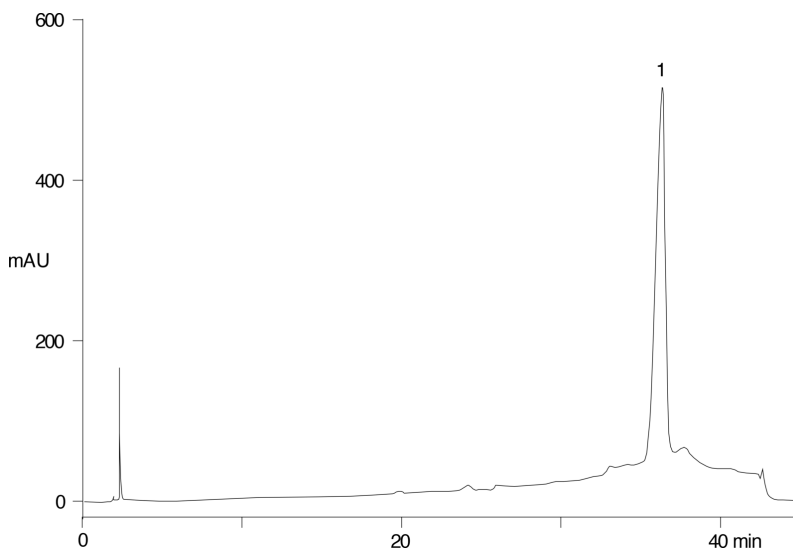

ANALYTES:

1 20mer Oligonucleotide



Products used in this application:

