

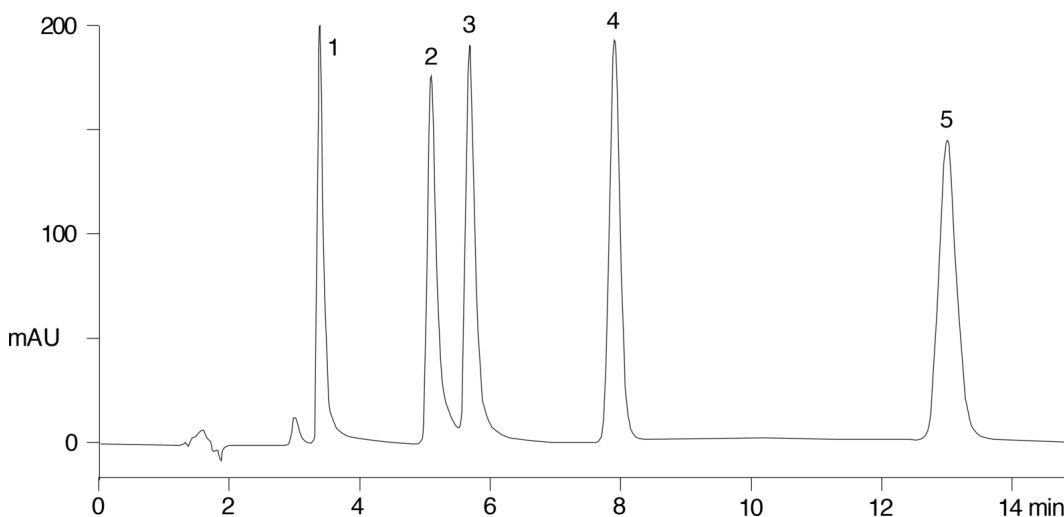
Benzodiazepams on Hyperclone 5u C18**Column:** HyperClone™ 5 µm ODS (C18) 120 Å, LC Column 150 x 4.6 mm, Ea**Dimensions:** 150 x 4.6 mm ID**Order No:** 00F-4361-E0**Elution Type:** Isocratic**Eluent A:** 20mM KH₂PO₄, pH=3.00, 400ul/L Triethylamine**Eluent B:** CH₃CN

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	63	37

Flow Rate: 1 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 230 nm (ambient)**Analyst Note:** SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, C18 4 x 3.0mm, 10/Pk Part No.: AJ0-4287

- Holder Part No.: KJ0-4282

**Products used in this application:****HyperClone™****ANALYTES:**

- 1 Bromazepam
- 2 Lorazepam
- 3 Oxazepam
- 4 Flunitrazepam
- 5 Diazepam

