

Nitro Aromatics on Hyperclone 5u Silica

Column: HyperClone™ 5 µm Silica 120 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4358-E0

Elution Type: Isocratic

Eluent A: Hexane

Eluent B: Ethyl Acetate

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	98	2

Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.


- SecurityGuard Cartridges, Silica 4 x 3.0mm, 10/Pk Part No.: AJ0-4348

- Holder Part No.: KJ0-4282



Products used in this application:

HyperClone™

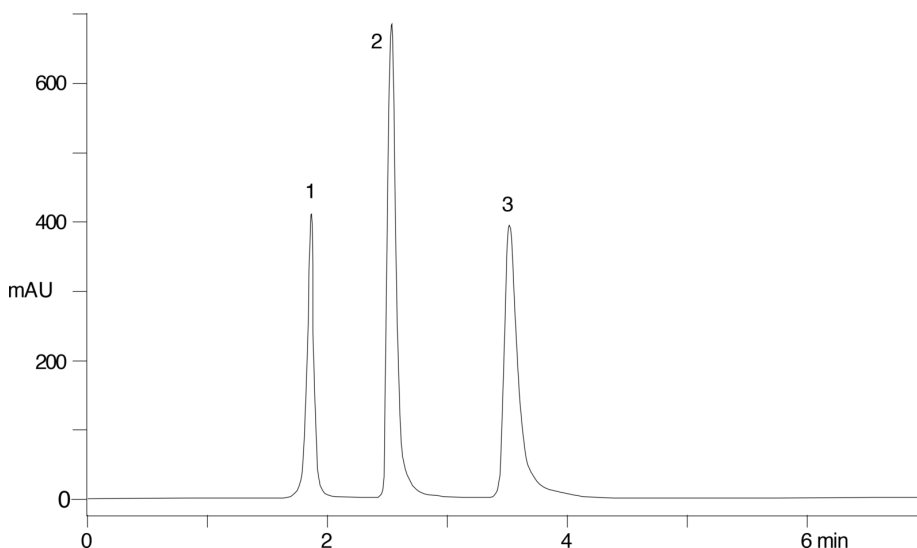

ANALYTES:

- 1 Toluene
- 2 Nitrobenzene
- 3 2,6-Dinitrotoluene

