

Proteins on Jupiter C4 300A

Column: Jupiter® 5 µm C4 300 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4167-E0

Elution Type: Gradient

Eluent A: 0.1% TFA : H₂O

Eluent B: 0.08% TFA : ACN

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0.5	95	5
	2	20	20	80

Flow Rate: 1 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ nm (ambient)


Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, Widespore C4 4 x 3.0mm, 10/Pk Part No.: AJ0-4330

- Holder Part No.: KJ0-4282

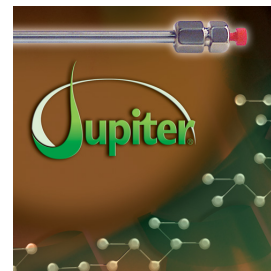

Products used in this application:

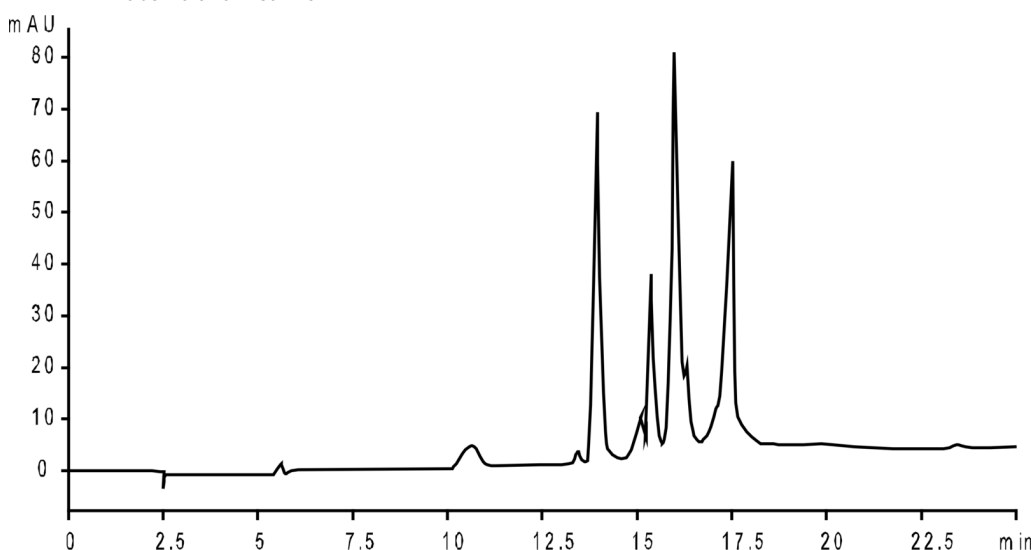

ANALYTES:

- 1 BSA
- 2 Glutamic dehydrogenase
- 3 B-Galactosidase
- 4 Ovalbumin

