

Tryptic Digest of Cytochrome c on Luna SCX

Column: Luna® 5 µm SCX 100 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4398-E0

Elution Type: Gradient

Eluent A: 20mM KPhos pH 2.5 + 25% Acetonitrile

Eluent B: A + 350mM KCl

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	1	100	0
	2	50	0	100

Flow Rate: 1 mL/min

Col. Temp.: 30 °C

Detection: UV-Vis Abs.-Variable Wave.(UV) @ nm (ambient)


Analyst Note: SecurityGuard™ Guard Cartridge System extends column lifetime.

- SecurityGuard Cartridges, SCX 4 x 3.0mm, 10/Pk Part No.: AJ0-4308

- Holder Part No.: KJ0-4282



Products used in this application:

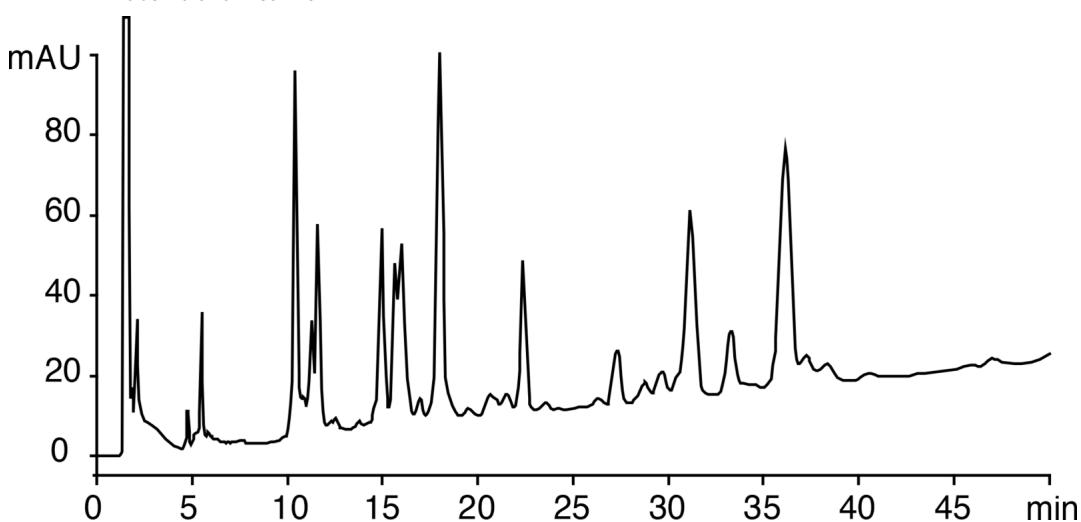

ANALYTES:

1 Cytochrome c

