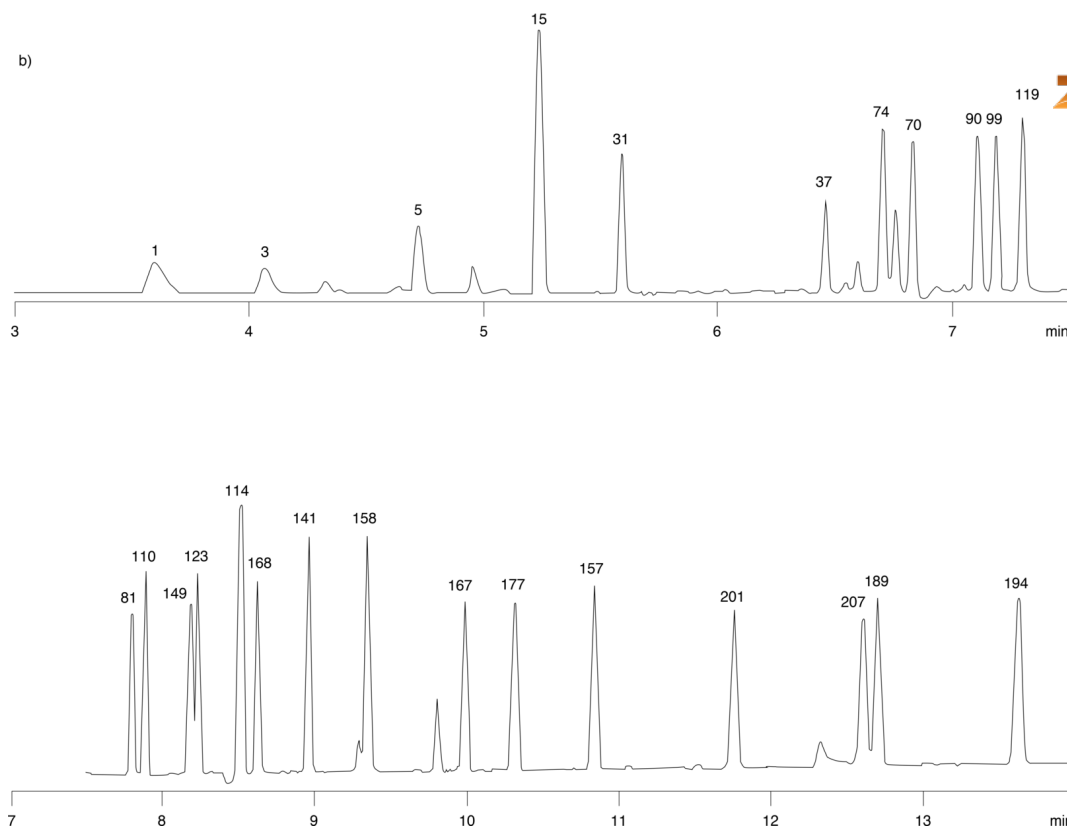


PCB Congener Analysis on ZB-1701 (Mix 2)

Column: Zebtron[™] ZB-1701, GC Cap. Column 30 m x 0.32 mm x 0.25 μ m, Ea
Phase: 14% Cyanopropylphenyl 86% Dimethylpolysiloxane
Dimensions: 30 meters x 0.32 mm x 0.25 μ m
Order No: 7HM-G006-11
Oven Profile: 130°C to 230°C at 20°C/min to 270°C at 4°C/min to 300°C at 20°C/min (hold 1min)
Carrier Gas: Constant Flow Helium, 2.5 mL/min
Injection: splitless :1 1 μ L @ 225°C
Detection: Electron Capture (ECD) (325°C)

14776

**Products used in this application:**

PCB Congener Analysis on ZB-1701 (Mix 2)

ANALYTES:

1	1
2	3
3	5
4	15
5	31
6	37
7	74
8	70
9	90
10	99
11	119
12	81
13	110
14	149
15	123
16	114
17	168
18	141
19	158
20	167
21	177
22	157
23	201
24	207
25	189
26	194

