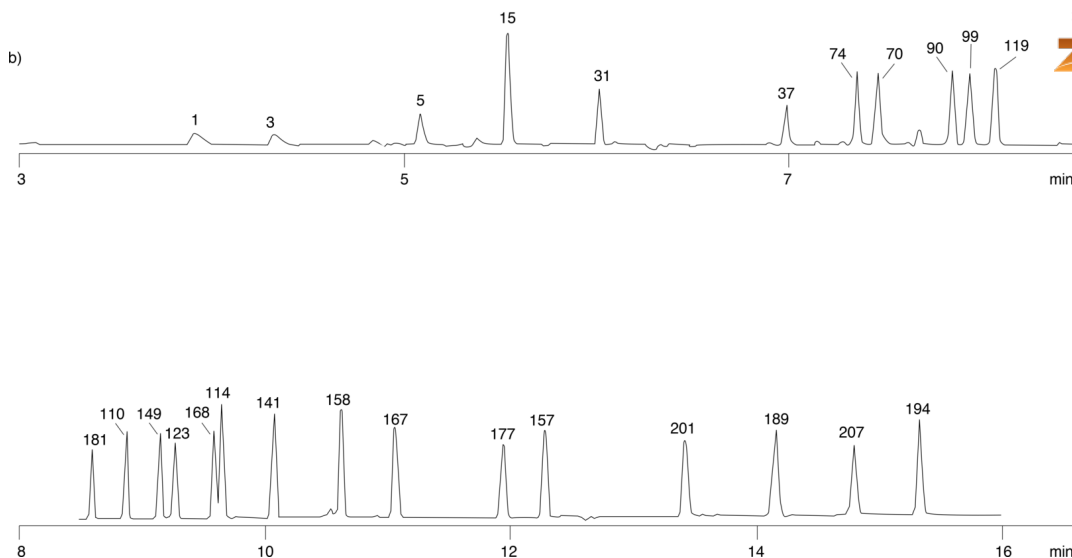


PCB Congener Analysis on ZB-50 (Mix 2)

Column: Zebtron[™] ZB-50, GC Cap. Column 30 m x 0.32 mm x 0.25 μ m, Ea
Phase: 50% Phenyl 50% Dimethylpolysiloxane
Dimensions: 30 meters x 0.32 mm x 0.25 μ m
Order No: 7HM-G004-11
Oven Profile: 130 °C to 230°C at 20 °C/min to 270 °C at 4 °C/min to 300 °C at 20 °C/min (hold 1 min)
Carrier Gas: Constant Flow Helium, 2.5 mL/min
Injection: Splitless :1 1 μ L @ 225°C
Detection: Electron Capture (ECD) (325°C)

14777

**Products used in this application:**

PCB Congener Analysis on ZB-50 (Mix 2)

ANALYTES:

1	1
2	3
3	5
4	15
5	31
6	37
7	74
8	70
9	90
10	99
11	119
12	181
13	110
14	149
15	123
16	168
17	114
18	141
19	158
20	167
21	177
22	157
23	201
24	189
25	207
26	194

