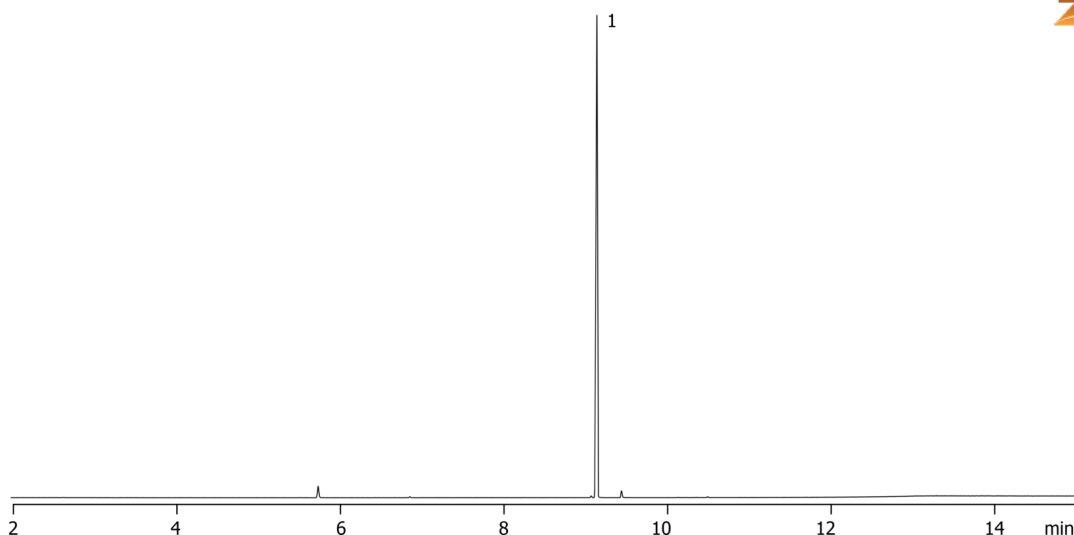


**PCP on ZB-5****Column:** Zebron<sup>™</sup> ZB-5, GC Cap. Column 15 m x 0.25 mm x 0.25 µm, Ea**Phase:** 5% Phenyl 95% Dimethylpolysiloxane**Dimensions:** 15 meters x 0.25 mm x 0.25 µm**Order No:** 7EG-G002-11**Oven Profile:** 70°C to 120°C at 10°C/min to 280°C at 20°C/min for 2 min**Carrier Gas:** Constant Flow Helium, 1.1 mL/min**Injection:** Split 10:1 1 µL @ 250°C**Detection:** Mass Selective (MSD) (180°C)**Analyst Note:** PCP at 133ng on column.

14786

**Products used in this application:****ANALYTES:****1** PCP