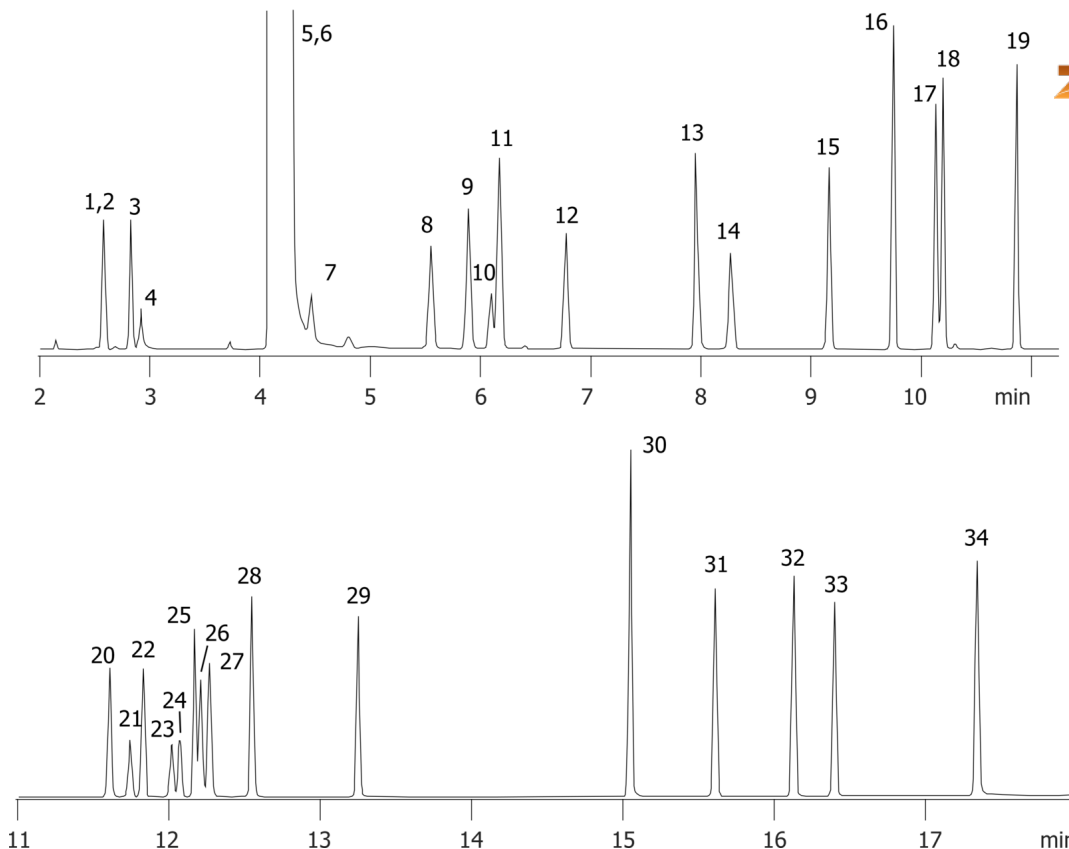


Alcohols on ZB-WAX 30x.32x0.5**Column:** Zebron[™] ZB-WAX, GC Cap. Column 30 m x 0.32 mm x 0.50 μ m, Ea**Phase:** Polyethylene Glycol (PEG)**Dimensions:** 30 meters x 0.32 mm x 0.5 μ m**Order No:** 7HM-G007-17**Oven Profile:** 60°C for 4 min to 90°C at 7°C/min to 220°C at 15°C/min.**Carrier Gas:** Constant Flow Helium, 1.2 mL/min**Injection:** Split 20:1 1 μ L @ 220°C**Detection:** Mass Selective (MSD) (180°C)

14795

**Products used in this application:**

Alcohols on ZB-WAX 30x.32x0.5

ANALYTES:

- 1 Methanol
- 2 tert-Butanol
- 3 Isopropanol
- 4 Ethanol
- 5 Chloroform
- 6 sec-Butanol
- 7 Propanol
- 8 2-Methylpropanol
- 9 3-Pentanol
- 10 Allyl alcohol
- 11 2-Pentanol
- 12 Butanol
- 13 3-Hexanol
- 14 Methallyl alcohol
- 15 Pentanol
- 16 4-Heptanol
- 17 Cyclopentanol
- 18 2-Ethyl-1-butanol
- 19 Hexanol
- 20 Cyclohexanol
- 21 cis-2-Methylcyclohexanol
- 22 trans-2-Methylcyclohexanol
- 23 cis-3-Methylcyclohexanol
- 24 cis-4-Methylcyclohexanol
- 25 Heptanol
- 26 trans-3-Methylcyclohexanol
- 27 trans-4-Methylcyclohexanol
- 28 2-Ethyl-1-hexanol
- 29 Octanol
- 30 Decanol
- 31 Phenylethanol
- 32 Phenylmethanol
- 33 2-Phenylethanol
- 34 3-Phenyl-1-propanol

