


Blood Alcohol on ZB-624**Column:** ZebronTM ZB-624, GC Cap. Column 30 m x 0.53 mm x 3.00 μ m, Ea**Phase:** 6% Cyanopropylphenyl 94% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.53 mm x 3 μ m**Order No:** 7HK-G005-36**Oven Profile:** 40°C for 6 min.**Carrier Gas:** Constant Flow Helium, 6.4 mL/min**Injection:** Headspace :1 1 μ L @ 225°C**Detection:** Flame Ionization (FID) (325°C)**Analyst Note:** Headspace analysis**Products used in this application:**

14806

**ANALYTES:**

- 1 Acetaldehyde
- 2 Methanol
- 3 Ethanol
- 4 Acetone
- 5 Isopropanol
- 6 n-Propanol

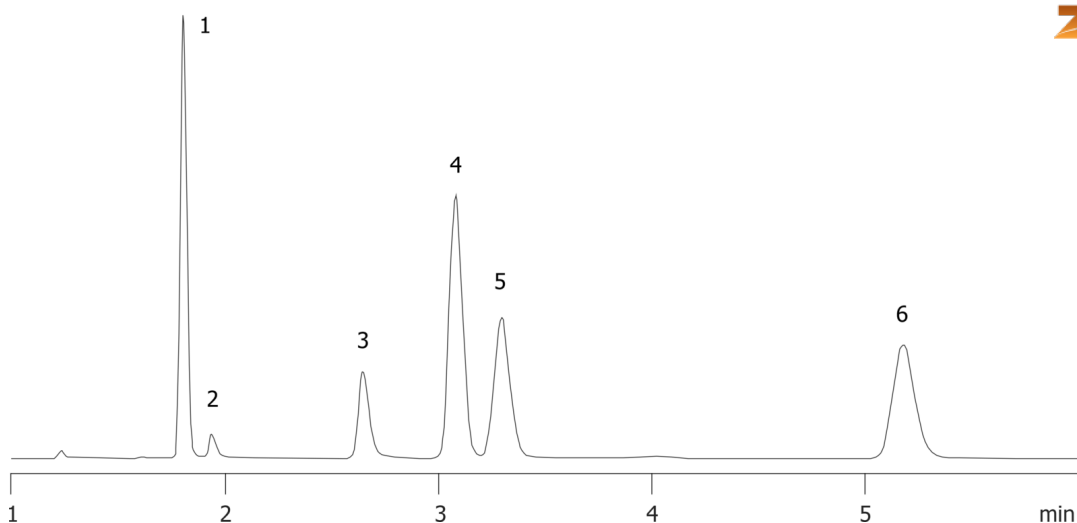