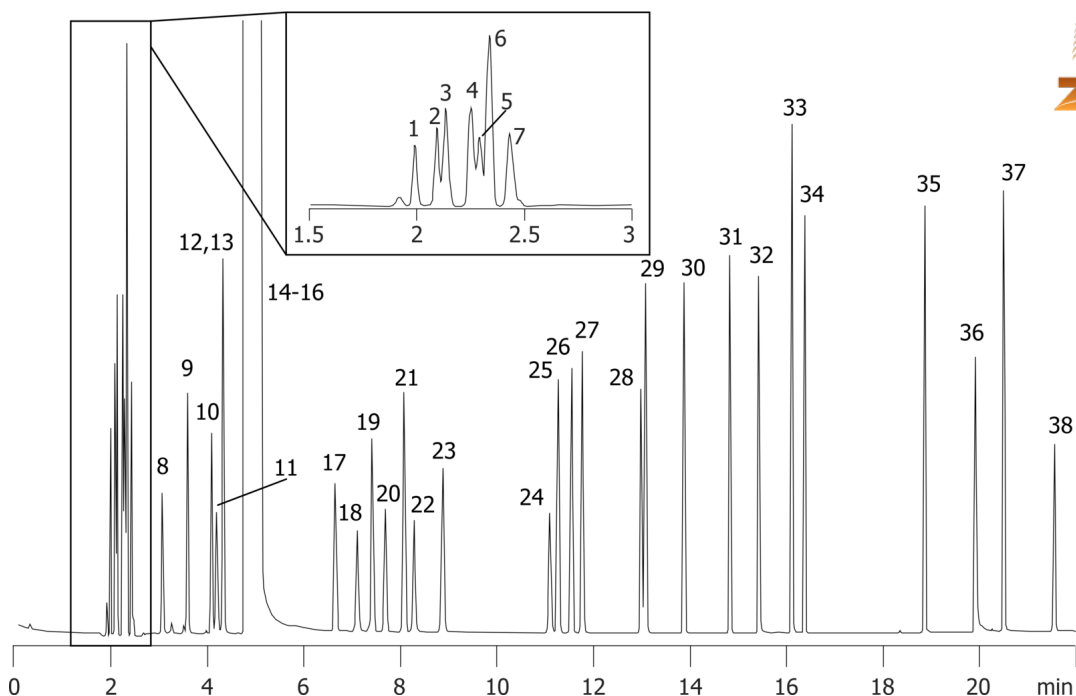


Industrial Solvents on ZB-WAX 30x0.25x0.25

Column: Zebtron™ ZB-WAX, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea
Phase: Polyethylene Glycol (PEG)
Dimensions: 30 meters x 0.25 mm x 0.25 µm
Order No: 7HG-G007-11
Oven Profile: 35°C for 5 min to 60°C at 4°C/min to 180°C at 12°C/min for 4 min.

Carrier Gas: Constant Flow Helium, 0.7 mL/min
Injection: split 100:1 1 µL @ 250°C
Detection: Mass Selective (MSD) (180°C)

14808

**Products used in this application:**

Industrial Solvents on ZB-WAX 30x0.25x0.25**ANALYTES:**

- 1** Pentane
- 2** Hexane
- 3** Ethyl ether
- 4** Isooctane
- 5** Methyl tert-butyl ether
- 6** Heptane
- 7** Cyclohexane
- 8** Acetone
- 9** Tetrahydrofuran
- 10** Ethyl acetate
- 11** Methanol
- 12** t-Butanol
- 13** MEK
- 14** Methylene chloride
- 15** Ethanol
- 16** Benzene
- 17** Trichloroethylene
- 18** Acetonitrile
- 19** Tetrachloroethylene
- 20** Chloroform
- 21** Toluene
- 22** 1-Propanol
- 23** 1,4-Dioxane
- 24** Allyl alcohol
- 25** Ethylbenzene
- 26** p-Xylene
- 27** m-Xylene
- 28** Pyridine
- 29** o-Xylene
- 30** Chlorobenzene
- 31** Styrene
- 32** 3-Picoline
- 33** N,N-Dimethylformamide
- 34** Anisole
- 35** Benzaldehyde
- 36** Ethylene glycol
- 37** N-Methyl-2-pyrrolidinone
- 38** Formamide

